

LM Dental™



LM

feel the difference



PRODUCT CATALOG

LM-DENTAL™

You can feel the difference

We are LM-Dental™, a dental instrument manufacturer located in Parainen, a town in southwest Finland.

Ever since founded in 1973 - more than 50 years ago - we have focused on developing, manufacturing, and marketing high-tech dental hand instruments. We apply the same high standards also to all other LM™ product groups: dental instrument tracking system, ultrasonic devices, and orthodontic appliances.

Our goal at LM™ is to create genuine value for dental care professionals by making instruments that are safe, efficient, and pleasant to use. We aim to make a difference where it counts: by helping people in all aspects of dentistry.

At LM™ we are truly dedicated to manufacturing excellence, allowing us to serve our customers with products that are precise, durable, and reliable. All our products are manufactured in Finland.

Continuous improvement in services and standards has always been our guiding principle, and by following that we've grown from a small garage workshop in the Finnish archipelago to a distinguished player in the global dental industry. Today we are the market leader in Nordic countries exporting 80% of our products globally.

WITH US YOU CAN FEEL THE DIFFERENCE.

Read more at www.lm-dental.com/About us/Who we are



PLANMECA

LM-Dental™ is part of the Planmeca Group,
one of the world's leading manufacturers of dental equipment.
www.planmeca.com

2

LM-DENTAL™

Feel the difference | LM™ product key features | Instrument handles and tip coatings

15

DIAGNOSTICS

Mirror handles | Explorers | Periodontal probes | Furcation probes

22

PERIODONTICS

Sickle scalars | Dual Gracey instruments | Universal curettes | Gracey curettes | Sharpen-free LM Sharp Diamond™ coated instruments | Diamond files | Hoe scalars | LM-ErgoMix™ implant instruments

41

RESTORATION

Cavity preparation instruments | Instruments for contouring and sculpting | Non-stick LM Dark Diamond™ coated instruments | LM-Arte™ instruments for esthetic restorations | LM-MultiHolder™ system

54

ENDODONTICS

Endodontic hand instruments

55

EXTRACTION

Elevators | Luxating instruments | Root tip instruments

58

ORTHODONTICS

Orthodontic hand instruments

59

INSTRUMENT HANDLING & MAINTENANCE

LM Dental Tracking System™ | LM-Servo™ instrument cassettes | LM-ServoMax™ tray system | Coding products | LM-RondoPlus™ sharpening machine

66

ORTHODONTIC APPLIANCES

LM-Activator™ | LM-Activator™ 2 | LM-OrthoSizer™ | LM-Trainer™



Making dental care safe, efficient and pleasant.

PIONEERING SINCE 1973

Instrument ergonomics

HOW TO DEFINE AN ERGONOMIC HAND INSTRUMENT?

Firstly a non-slip, natural grip on an instrument are essential. Secondly a lightweight and optimized design prevents effectively occupational diseases. Thirdly a good tactile sensitivity and easy handling of the instrument guarantee the best clinical results. Reliability, wear resistance and comfort are the features that will be evident for both, the dentist and the patient.

- These are the guiding principals how Pekka Kangasniemi, DDS, one of the founders of LM-Dental™, defined instrument ergonomics already in the 1980s.

PIONEERING INNOVATIONS IN DENTAL INSTRUMENT ERGONOMICS

As ergonomics being one of the first key concepts for LM™ product line, were the oversized handles the first step on that path and introduced already in 1974. In the end of 1980s ergonomics were further enhanced by the first silicone handles and the year 1995 marked the year for our first oversized silicone handle LM-ErgoMax™.

The greatest outcome of intense research and cooperation with dental clinicians was the launch of LM-ErgoSense® instruments in 2015 with even more ergonomic grip and with possibility to integrate a trackable RFID tag into the handle making it the first "smart instrument" on the market.

What LM™ brand is widely known and recognized is the innovation introduced already in 1980: the color-coded handles to ease the identification of the instrument. These radical product theses at the time - a thicker handle, optimized shape, elastic surface, lightness and color-coding - have become industry wide standards in the 21st century. We are the innovation leader that others copy.



50 YEARS OF MAKING
THE DIFFERENCE

Read more at
www.lm-dental.com





Fresh color-coding

- ease the identification during clinical procedures and maintenance

Elastic, ergonomic silicone surface and optimal contouring

- comfortable and non-slip grip

Optimal rigidity

- excellent sensitivity and efficiency
- outstanding corrosion resistance

Special alloyed metal core

- lightweight structure
- metal-to-metal connection for enhanced tactile sensitivity

LM-ErgoSense® handle *)

- oversized large diameter
- available with integrated LM DTS™ RFID tag
- available with sharpen free or non-stick coated tips

Precision formed tips and high-quality finishing

- extremely sharp and smooth cutting edge

LM-DuraGradeMax™ stainless steel tips

- high wear resistance
- immune to corrosion
- uncoated: can be resharpened

*) See the other handle options on pages 12-13.



MAKING DENTAL CARE SAFE, EFFICIENT AND PLEASANT.

We are truly dedicated to manufacturing excellence, allowing us to serve our customers with products that are precise, durable, and reliable. Continuous improvement in services and standards has always been our guiding principle. We apply the same high standards to our all product categories and are number-one-partner with the best-in-class dental products and smart solutions.

LM-ErgoSense® instrument handle

ENHANCING COMFORT AND PERFORMANCE

When holding the instruments for approximately 2,000 hours each year, the design, weight and feel of the instrument are crucial for the dental practitioner in terms of his/her performance and productivity. An ergonomic, high quality hand instrument naturally feels good in the hand.

The LM-ErgoSense® instrument handle is outcome of intense research and product development in close cooperation with dental clinical experts. It is truly a unique combination of ergonomic, functional design and high-tech in form of integrated RFID technology.

OVERSIZED Ø 13.7 MM SILICONE HANDLE

A hand instrument with thick silicone handle has been found to be more usable, cause lower perceived strain and to be more productive than those with thinner handles (1).

RATED AS BEST HAND INSTRUMENT


LM-ErgoSense® was rated as the best in 10 out of 14 specified usability features compared with the other tested instruments in a dental scaling simulation (1).

"FEELS GREAT IN THE HAND"

In study focusing on ergonomics, usability characteristics and work productivity of different hand instruments the thick, silicone coated LM™ instruments were the most preferred choice in regard to all 18 usability criteria used in the test (2).

UNIQUE INTELLIGENCE INSIDE

Integrated LM DTS™ RFID tag offers unique opportunities to improve material handling and patient safety. Read more on page 60.

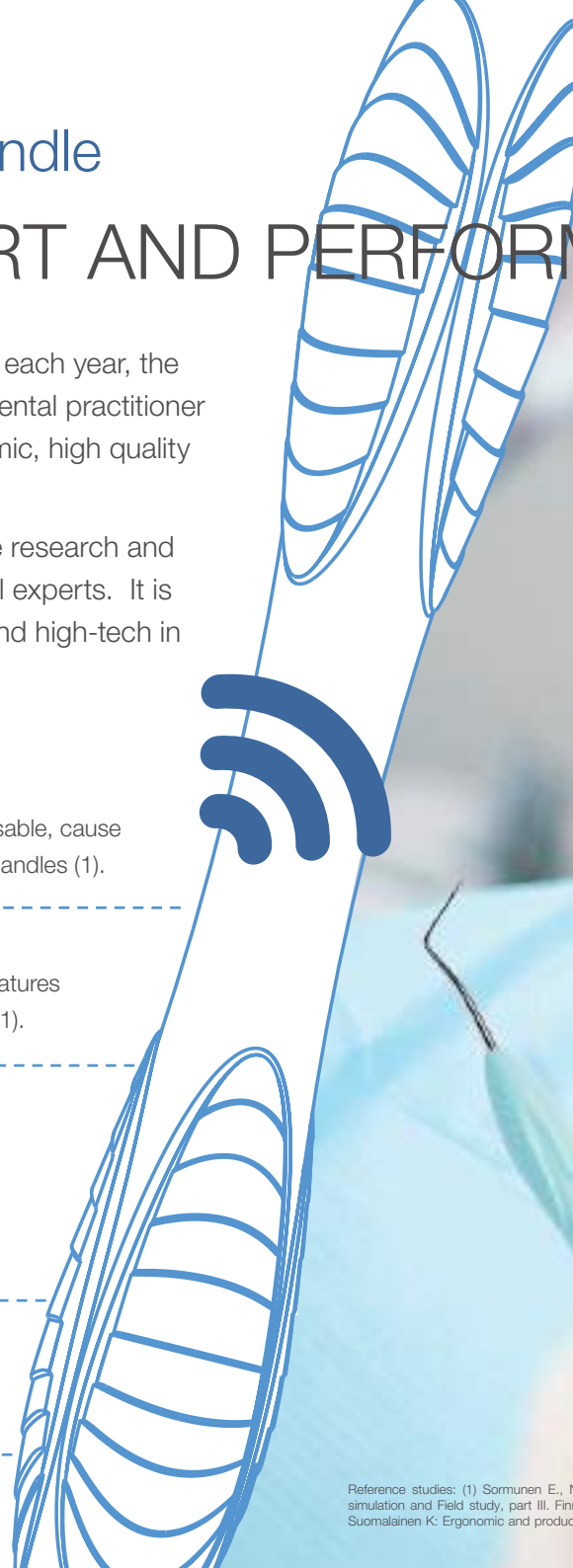

LM SHARP DIAMOND™


Sharpen-free instruments

Sharpen free LM Sharp Diamond™ coated instruments are ideal choice for clinicians who frequently use hand instruments and want to minimize the time spent on instrument maintenance. These instruments present optimal characteristics of superior sharpness, tactile sensitivity and comfort while combined with premium ergonomics of LM-ErgoSense® handle.

LM Sharp Diamond™ instrument tip is made from exceptionally durable special metal alloy, and its wear resistance is enhanced by a protective micro membrane coating providing excellent sharpness and control. LM Sharp Diamond™ instruments are well-suited for all forms of scaling and tips retain their original functionality, size and shape throughout their lifetime. The time saved from sharpening can be used for productive work and patient care.

What do the dental hygiene experts say? Scan the QR for the video.

See the periodontal product range on pages 22-40.

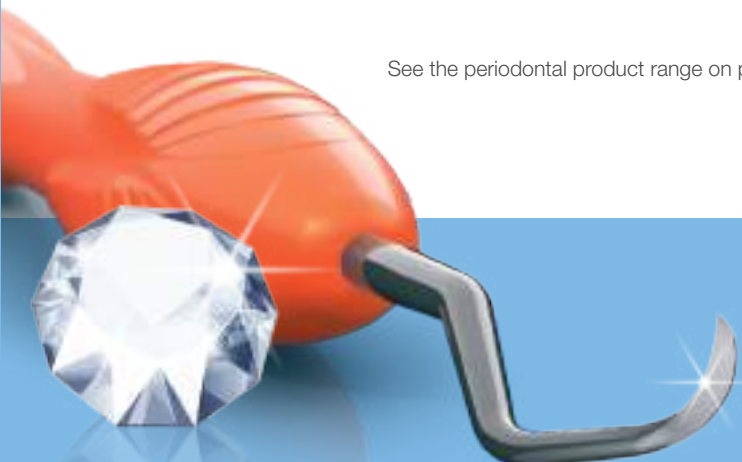

With LM Sharp Diamond™, I'm not just doing my job - I'm perfecting it. I choose to work smart. I choose to work efficiently. I choose LM Sharp Diamond™."


Dr. Joy Void-Holmes

RHD, BSDH, DHSc, CEO and founder of Dr. Joy




CHOOSE SHARP.

SHARPEN FREE MICRO COATING

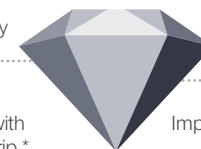Savings of time & money Long-lasting superior sharpness
 Highest rated ergonomics with sensational grip * Improved productivity & clinical results



DIAMOND-LIKE FEATURES

Sharpen free 
 Superior sharpness 
 Highest rated ergonomics..... 
 Minimum maintenance..... 


HARDNESS OVER
63
 HRC

HARDNESS

Hardness contributes to the durability of the instrument. LM Sharp Diamond™ reaches hardness of over 63 on Rockwell scale.

WEAR RESISTANCE **

82% decrease in wear compared to typical dental hand instrument material.

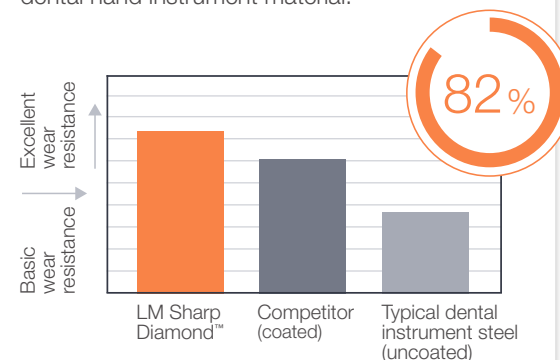




About SHARPNESS of dental hygiene instruments

The typical lifespan of a dental instrument varies between 6 to 18 months, depending on usage frequency and handling techniques. For optimal performance, most hygienists recommend sharpening instruments at least every two weeks. Prolonged usage without regular sharpening can result in bluntness, reducing the effectiveness of your work and ultimately diminishing patient comfort and satisfaction.

A significant number of instruments are often used beyond the recommended time frame, leading to decreased performance due to dullness. While sharpening is vital for maintaining sharpness, excessive sharpening can alter the original cutting shape and angularity, potentially decreasing instrument durability and increasing the risk of breakage.

It is of utmost importance for hygienists to utilize sharp instruments. Sharp scalers not only reduce fatigue and ergonomic strain for hygienists but also contribute to enhanced patient comfort.

COMPARISON of sharpen free instruments and instruments to be sharpened

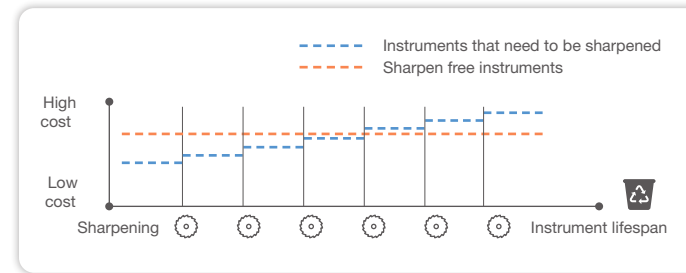


to easily 5-10 € per sharpening. As a result, the price difference between sharpen free hand instruments and those that require sharpening becomes negligible already after approximately 3-4 sharpening sessions.

By choosing sharpen free hand instruments, the need for sharpening is eliminated. This not only saves valuable time but also reduces expenses, which can then be redirected towards more productive work.

COST PERSPECTIVE

Switching from conventional to sharpen free instruments a small clinic allows to reduce costs significantly, even 10 000 €, in one year.


The purchase price of sharpen free instruments is ~25% higher compared to instruments that can be reshaped. However, the cost of sharpening services or doing it yourself, counting time, amounts

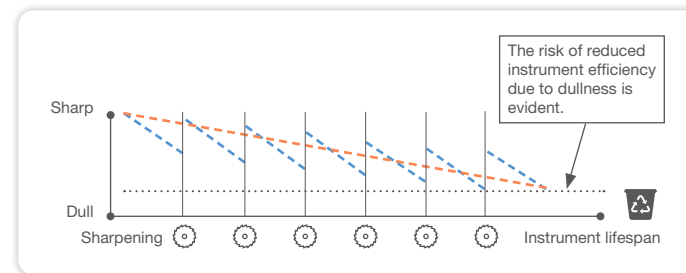

SHARPNESS AND EFFECTIVENESS PERSPECTIVE

Sharpen free instruments can retain their sharpness for 82% longer compared to instruments that are reshaped, which tend to lose their sharpness more quickly.

Additionally, reshaping instruments display greater fluctuations in sharpness over time.

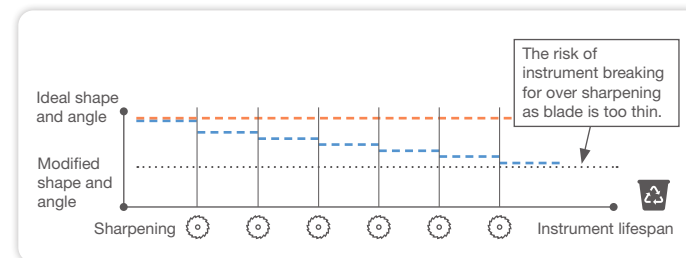

SHAPE AND SAFETY PERSPECTIVE

Sharpen free instruments maintain their original design, shape, and angulation consistently throughout their entire lifespan.

Excessive sharpening of instruments can alter the original cutting shape and angularity, which may lead to reduced instrument durability and increased risk of breakage.

LM DARK DIAMOND™

Non-stick coated instruments



LM Dark Diamond™ instruments with non-stick coating are designed for easy placing and contouring of composite resins. The diamond-like coating reduces the stickiness and allows pleasant and precise restorative instrumentation.

- Contemporary diamond-like carbon (DLC) coating reduces stickiness for smooth composite placement
- Hard surface protecting the tip from scratches
- Smooth surface for easy cleaning
- Dark coating color for minimal reflections, and maximal contrast against filling
- Ergonomic LM-ErgoSense® handle for excellent precision and feel

LM-Arte™ Dark Diamond edition

PLEASANT AND PRECISE SCULPTING OF COMPOSITE

Non-stick coated LM-Arte™ Dark Diamond instruments are the latest outcome of StyleItaliano™ and LM-Dental™ collaboration featuring exclusive coating for smooth and precise non-stick placement and sculpting of composite resins.

What do the StyleItaliano™ experts say?
Scan the QR for the video.

See LM-Arte™ product range on pages 45-48.



The sensitivity, flexibility and precision are the same as with conventional LM-Arte™ instruments, but the LM Dark Diamond™ coating eliminates the issue of resin sticking to spatulas, even for the most delicate resins.”

Dr. Paulo Monteiro Honorary Member of StyleItaliano™



LM-Dental™ and Styleitaliano™

LM-Arte™ restorative instruments

SIMPLIFYING ESTHETIC RESTORATIONS

LM-Arte™ is a range of innovative instruments for esthetic restorations that is designed especially for composite layering. Each instrument is color-coded and named after its main function, enabling you to administer easy, efficient and quick restorative treatments. The lightweight and ergonomic LM™ handle design, together with high-quality tip materials, provides excellent tactile sensitivity and easy handling or the best clinical results.

CO-OPERATION MAKING THE DIFFERENCE

LM-Arte™ instrument series has been developed in collaboration with Styleitaliano™, a group of skilled dentists who share our passion for improving dentistry.

The Styleitaliano™ movement emerged from a simple idea of two friends - Dr. Walter Devoto and Dr. Angelo Putignano. The renowned experts in conservative and esthetic dentistry had noticed that although materials suitable for creating beautiful direct restorations were available, appropriate instruments for employing these successfully were lacking. With LM-Dental™'s dedication to manufacturing excellence the unique concept of LM-Arte™ instruments was developed, and it has been a breakthrough in esthetic dentistry.



BE PART OF THE SUCCESS!

Visit Styleitaliano™ website at www.styleitaliano.org for restorative clinical cases with LM-Arte™ instruments.



Follow and join the social media channels!



Choosing the handle

Choosing the most ergonomic hand instrument can make a big difference in terms of discomforts and disorders when performing dental work. It is by no means insignificant what is the diameter or material of the grip surface of the hand instrument handle.

LM™ hand instruments are available with a selection of different handles all of them exceptionally ergonomic due to the elastic silicone coating and optimized contouring. The clinician's hand fit and indication both play a role in selecting an ideal handle type. Interchanging different handles, may place the user's hand under even less strain.



COLORS, NAMING AND CODING

To facilitate color coding and brighten up working practices all the LM™ instruments come in refreshing colors. The name and the code number on the handle ease the identification of the instrument during clinical procedures and maintenance.

Handle selection



LM-ErgoSense®

The larger diameter and design provide sensational grip that has been proven to be ergonomic and efficient in scientific clinical testing.

LM-ErgoSense® handle is available with integrated LM DTS™ RFID tag. The compatibility with the LM Dental Tracking System™ offers many opportunities to improve asset management and patient safety.



LM-ErgoMax™

Classic ergonomic design that others copy. Functions well in all clinical procedures and guarantees excellent tactile sensitivity and relaxed grip.



Codes, symbols and diameters

LMErgoSense®

ES = LM-ErgoSense® handle
(Ø 13.7 mm)



EST = LM-ErgoSense® handle with integrated LM DTS™ RFID tag

LMErgoMax™

XSI = LM-ErgoMax™ handle
(Ø 11.5 mm)

Special handles



LM-ErgoNorm™

The thinner handle that functions best in light clinical procedures.



LM-ErgoMix™

The LM-ErgoMix™ handle for exchangeable tips offers all the economical and ecological benefits of retipping. See page 39.



LM-ErgoSingle™

Well balanced mirror handle with better reachability. See page 15.

LMErgoNorm™

SI = LM-ErgoNorm™ handle
(Ø 8.5 mm)

LMErgoMix™

EM = LM-ErgoMix™ handle
(Ø 12.5 mm)

LMErgoSing™

ESI = LM-ErgoSingle™ handle
(Ø 11.5 mm)

Choosing instruments with special coated tips

SHARPEN FREE PROPHY AND PERIO

Sharpen free LM Sharp Diamond™ instrument tips are made from exceptionally durable special metal alloy, and their wear resistance is enhanced by a protective micro membrane coating made using the most advanced PVD (Physical Vapor Deposition) technology. Choose sharp and work efficiently.



LMErgoSense®

SDES = LM-ErgoSense® (ES) handle with sharpen free LM Sharp Diamond™ (SD) coating



NON-STICK RESTORATIVE INSTRUMENTS

LM Dark Diamond™ instruments with non-stick coating allow easy placing and contouring of composite resins. Dark colored tips provide maximum contrast against the filling and reduce the reflections.



LMErgoSense®


DDES = LM-ErgoSense® (ES) handle with non-stick LM Dark Diamond™ (DD) coating

PATIENT SAFETY AND TRACEABILITY

Coated instruments, like all LM™ hand instruments with LM-ErgoSense® handle, are available with integrated LM DTS™ RFID tag for the compatibility with unique LM Dental Tracking System™.



SDEST / DDEST = LM-ErgoSense® (ES) handle with sharpen free (SD) or non-stick (DD) coating and integrated LM DTS™ RFID tag (T)

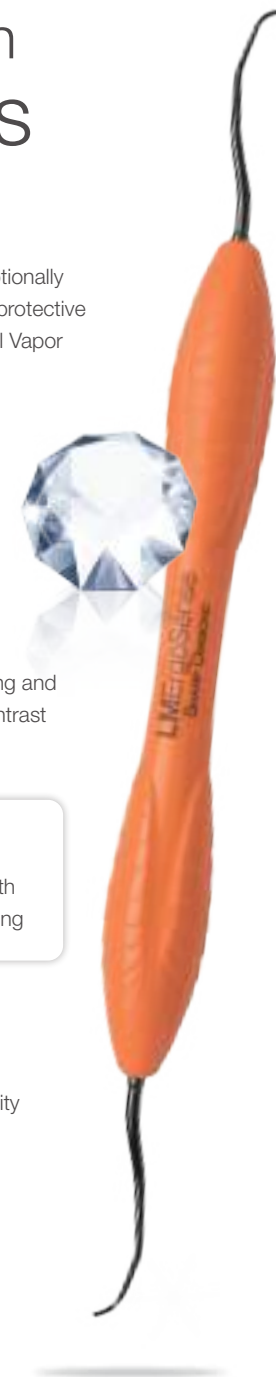



Making dental care safe, efficient and pleasant.

MIRROR HANDLES

Simple stem

LM 25 mirror handle threads are in accordance with international standard M2.5. The mirror shank must have a diameter of 2.5 mm.

Cone socket

LM 28 mirror handle threads are in accordance with international standard UNF4-48. The mirror shank must have a diameter of 2.6 mm.



LM-ErgoSingle™ Mirror Handle

- Well balanced, better reachability
- Color options: red, grey and blue

LMErgoSingle™

- 25ESIRED
- 25ESIGREY
- 25ESIBLUE

LMErgoSingle™

- 28ESIRED
- 28ESIGREY
- 28ESIBLUE



Mirror Handle 25

- Simple stem design, single-ended
- Color options: red, grey and blue

LMErgoSense®

- 25ESRED
- 25ESGREY
- 25ESBLUE

LMErgoSense®

- 25ESREDT
- 25ESGREYT
- 25ESBLUET

LMErgoMax™

- 25XSIREDD
- 25XSIGREYD
- 25XSIBLUED

LMErgoNorm™

- 25SIREDD
- 25SIGREYD
- 25SIBLUED

i
Handle options:

LMErgoSense® (Cone socket)

LMErgoMax™ (Simple stem)

LMErgoNorm™ (Simple stem)



Mirror Handle 28

- Cone socket design, UNF4-48 threads, single-ended
- Color options: red, grey and blue

LMErgoSense®

- 28ESRED
- 28ESGREY
- 28ESBLUE



- 28ESREDT
- 28ESGREYT
- 28ESBLUET

LMErgoMax™

- 28XSIRE
- 28XSIGREY
- 28XSIBLUE



Mirror Handle, ortho

- Orthodontic mirror handle with Bausch-Verbiest design
- For measuring the sagittal overjet and vertical overbite
- Millimeter scale, total length 25 mm

LMErgoSense® 25-26ES

LMErgoSense® 25-26EST

LMErgoMax™ 25-26XSI

LMErgoSense® 28-26ES

LMErgoSense® 28-26EST

LMErgoMax™ 28-26XSI



Mirror Handle JM

- For implant installation

LMErgoSense® 25-473ES

LMErgoSense® 25-473EST

LMErgoMax™ 25-473XSI



Handle options:



LMErgoSense®



LMErgoMax™



LMErgoNorm™

EXPLORERS

For diagnosing calculus and caries and exploring of pockets, restorations and furcations.



Explorer 5-8

- Standard explorer for diagnostic and operative procedures

	5-8ES
	5-8EST
	5-8XSI
	5-8SI




Explorer 5-9

- Standard explorer for diagnostic and operative procedures

	5-9ES
	5-9EST
	5-9XSI
	5-9SI

Explorer, centered 9-17C

- Standard explorer for diagnostic and operative procedures
- Especially for probing distal surfaces of posterior teeth

	9-17CES
	9-17CEST
	9-17CXSI
	9-17CSI



Explorer 17-23

- For distal surfaces (17)
- Elongated lower shank and "back-action" blade (17)
- For areas in the vicinity of proximal contact surfaces (17)
- Also for probing occlusal surface grooves (23)
- For detecting subgingival calculus (23)

	17-23ES
	17-23EST
	17-23XSI
	17-23SI

Explorer 17C-23





- For distal surfaces, centered tip (17C)
- For areas in the vicinity of proximal contact surfaces (17C)
- Also for probing occlusal surface grooves (23)
- For detecting subgingival calculus (23)

	17C-23ES
	17C-23EST
	17C-23XSI
	17C-23SI



Explorer 3-6

- Standard explorer for diagnostic and operative procedures

	3-6ES
	3-6EST
	3-6XSI
	3-6SI

Explorer 11-12

- Angles similar to Gracey 11/12 (211-212)
- For probing proximal and cervical caries
- For subgingival root examination

	11-12ES
	11-12EST
	11-12XSI
	11-12SI





i
Handle options:

LMErgoSense®

LMErgoMax™

LMErgoNorm™

Flexplorer 11F-12F

- Flexible explorer angulated similarly to Gracey 11/12
- Well-suited for periodontal diagnosing of subgingival calculus and probing deep pockets

	11F-12FES
	11F-12FEST
	11F-12FXSI

Explorer 18-19

- For prosthetic crowns, bridges and inlay-/onlay fillings
- Excellent sensitivity at seams

	18-19ES
	18-19EST
	18-19XSI
	18-19SI

Explorer 29

- Holst design
- Lateral angle on the shank
- Single-ended

	29ES
	29EST
	29XSI
	29SI

Explorer 9

- Single-ended

	9ES
	9EST
	9XSI
	9SI

PERIODONTAL PROBES

Combination instruments for probing caries and parodontium.

- Most periodontal probes are available also in single-ended designs.
- Color-coded scales facilitate measurement of the gingival pocket



2 mm scale




2 mm scale
ball end



2 mm scale
ball end



3 mm scale
ball end

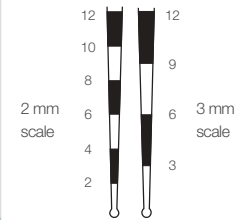

Ball end
ø 0,5 mm

- Increases tactile sensitivity
- Does not damage tissue or pierce pocket bottom
- Gives more accurate pocket depth measurement
- Gives information on subgingival calculus and other root surface unevenness

Explorer – Periodontal Probe 23-52B

- LM**ErgoSense® 23-52BES
- CC** 23-52BEST
- LM**ErgoMax™ 23-52BXS
- LM**ErgoNorm™ 23-52BS

Explorer – Periodontal Probe 8-520B


- LM**ErgoSense® 8-520BES
- CC** 8-520BEST
- LM**ErgoMax™ 8-520BXS
- LM**ErgoNorm™ 8-520BS

Explorer – Periodontal Probe 23-520B

- LM**ErgoSense® 23-520BES
- CC** 23-520BEST
- LM**ErgoMax™ 23-520BXS
- LM**ErgoNorm™ 23-520BS

Explorer – Periodontal Probe 23-530B

- LM**ErgoSense® 23-530BES
- CC** 23-530BEST
- LM**ErgoMax™ 23-530BXS
- LM**ErgoNorm™ 23-530BS

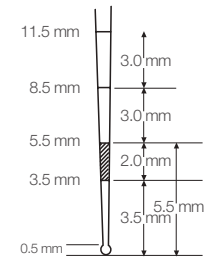

WHO probe was developed for CPITN method.

Diagnosis:

- 0 = healthy
- 1 = bleeding during probing
- 2 = supra- and/or subgingival calculus, filling or crown excesses
- 3 = gingival pockets to 5 mm
- 4 = pockets deeper than 6 mm

Williams Probe

- Single-ended

LMErgoSense®	51ES
☪	51EST
LMErgoMax™	51XSI
LMErgoNorm™	51SI

Explorer – WHO Probe 8-550B

LMErgoSense®	8-550BES
☪	8-550BEST
LMErgoMax™	8-550BXSI
LMErgoNorm™	8-550BSI


Explorer – WHO Probe 23-550B

LMErgoSense®	23-550BES
☪	23-550BEST
LMErgoMax™	23-550BXSI
LMErgoNorm™	23-550BSI

WHO Probe 550B

- Single-ended

LMErgoSense®	550BES
☪	550BEST
LMErgoMax™	550BXSI
LMErgoNorm™	550BSI

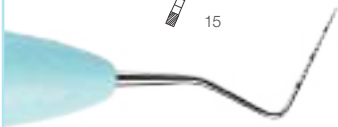

North Carolina Probe 54B

- Single-ended

LMErgoSense®	54BES
☞	54BEST
LMErgoMax™	54BXS



Marquis 12 Probe 53B

- Single-ended
- Similar to 530B but without ball end

LMErgoSense®	53BES
☞	53BEST
LMErgoMax™	53BXS

FURCATION PROBES



For probing the extent and depth of furcation lesions.



- For probing furcation lesions from different angles, both in the upper and in the lower jaw region.



Nabers Probe

LMErgoSense®	20-21ES
☞	20-21EST
LMErgoMax™	20-21XS
LMErgoNorm™	20-21SI

Furcation Probe

LMErgoSense®	20B-21BES
☞	20B-21BEST
LMErgoMax™	20B-21BXS
LMErgoNorm™	20B-21BSI



Handle options:



LM ErgoSense® ☞



LM ErgoMax™



LM ErgoNorm™



SICKLE SCALERS



For the removal of supragingival calculus.

Instrumentation:

For calculus removal place the tip third of the scaler's cutting edge against the tooth. Tilt the instrument toward the tooth to achieve 70° to 80° angle between the tooth and the blade. Apply lateral pressure and activate the scaler by using vertical, diagonal or horizontal pull strokes. To maintain control use short 2–3 mm long strokes.

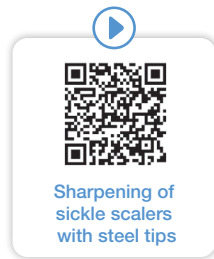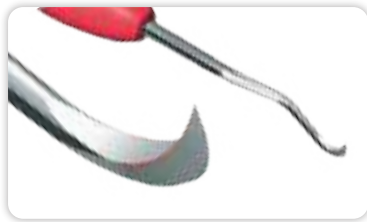

Design:

- Triangular cross section
- Face perpendicular to the lower shank
- Two cutting edges (red line), pointed tip



Handle options:



LM-SharpJack™ SD

- Anterior/posterior scaler
- Thin dual-edged blades for easy and efficient interproximal access
- With sharpen free LM Sharp Diamond™ coated tips



Sickle LM23 Slim SD

- Universal scaler for anterior and posterior areas
- Thin version of popular LM23 scaler (313-314)
- Better access to interproximal areas
- With sharpen free LM Sharp Diamond™ coated tips





Delicate working end. Especially for tight interdental spaces.

Micro Sickle

- Anterior/posterior scaler
- Delicate working end, elongated shank
- Ideal for tight interproximal surfaces

- LMErgoSense®** 301-302ES
- LMErgoSense®** 301-302EST
- LMErgoMax™** 301-302XSI
- LMErgoNorm™** 301-302SI



Micro Sickle SD

- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 301-302SDES
- LMErgoSense®** 301-302SDEST



For all tooth surfaces, especially proximal surfaces.

Mini Sickle

- Universal scaler
- For all tooth surfaces
- Fine, slightly angled blade

- LMErgoSense®** 311-312ES
- LMErgoSense®** 311-312EST
- LMErgoMax™** 311-312XSI
- LMErgoNorm™** 311-312SI



Mini Sickle SD

- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 311-312SDES
- LMErgoSense®** 311-312SDEST



Sickle LM23

- Posterior scaler
- Shorter blade and longer lower shank than in the Mini Sickle (311-312)

- LMErgoSense®** 313-314ES
- LMErgoSense®** 313-314EST
- LMErgoMax™** 313-314XSI
- LMErgoNorm™** 313-314SI



Sickle LM23 SD

- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 313-314SDES
- LMErgoSense®** 313-314SDEST





Handle options:



LM ErgoSense®

LM ErgoSense®



LM ErgoMax™



LM ErgoNorm™



Sickle NV1 SD

- Anterior scaler
- Sickle end for removing interproximal calculus and stain
- Disc end for lingual surfaces of incisors
- With sharpen free LM Sharp Diamond™ coated tips


LM ErgoSense® 317-318SDES
 LM ErgoSense® 317-318SDEST



Sickle NV2 SD

- Posterior scaler
- Heavily angled long shank allows for posterior access
- Thin and long blades provide easy access to interproximals
- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense® 319-320SDES
 LM ErgoSense® 319-320SDEST

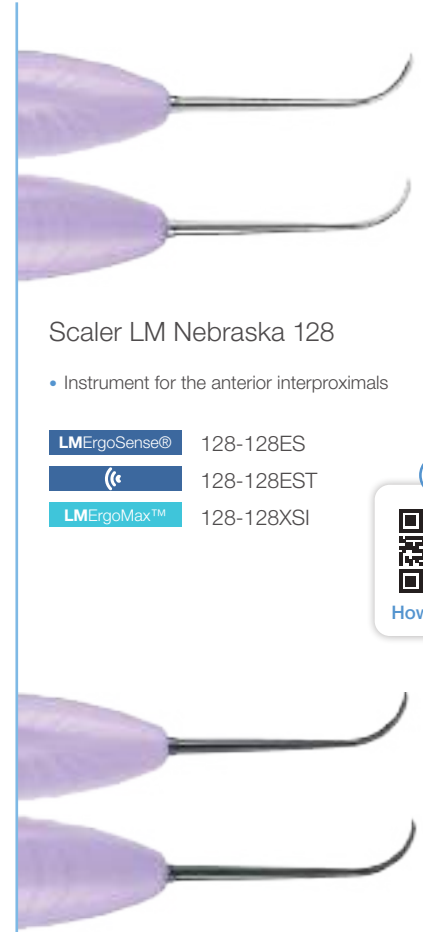

Scaler LM Nebraska 128


- Instrument for the anterior interproximals

LM ErgoSense® 128-128ES
 LM ErgoSense® 128-128EST
 LM ErgoMax™ 128-128XSI

Scaler LM Nebraska 128 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense® 128-128SDES
 LM ErgoSense® 128-128SDEST

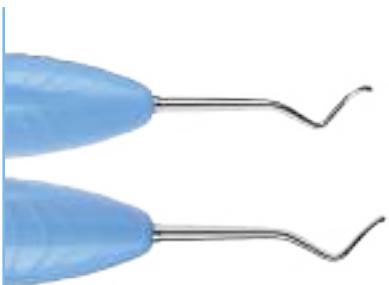

Sickle LM204S

- Sharply angulated shank and very fine working end

- LMErgoSense®** 303-304ES
- LMErgoSense®** 303-304EST
- LMErgoMax™** 303-304XSI
- LMErgoNorm™** 303-304SI



Scaler H6-H7

- For all tooth surfaces
- Longer shank and smaller lateral bending than in LM-SharpJack™ (315-316)

- LMErgoSense®** 146-147ES
- LMErgoSense®** 146-147EST
- LMErgoMax™** 146-147XSI
- LMErgoNorm™** 146-147SI



Sickle LM204S SD

- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 303-304SDES
- LMErgoSense®** 303-304SDEST



Scaler H6-H7 SD

- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 146-147SDES
- LMErgoSense®** 146-147SDEST



Push Scaler H4-H5

- The chisel-like working end used by pushing
- The jacquette working end used by pulling
- 90° working surface

- LMErgoSense®** 144-145ES
- LMErgoSense®** 144-145EST
- LMErgoMax™** 144-145XSI
- LMErgoNorm™** 144-145SI



Laplander

- Similar to H6-H7 (146-147) but with a lateral bending
- Easy access to molar area and proximal surfaces

- LMErgoSense®** 148-149ES
- LMErgoSense®** 148-149EST
- LMErgoMax™** 148-149XSI



Crane-Kaplan

- Extremely strong, rigid, sharply angled blade
- For removal of heavy supragingival calculus

- LMErgoSense®** 322-323ES
- LMErgoSense®** 322-323EST
- LMErgoMax™** 322-323XSI
- LMErgoNorm™** 322-323SI



Mini Kaplan

- Same angle but with a more delicate blade than in the Crane-Kaplan (322-323)

- LMErgoSense®** 333-334ES
- LMErgoSense®** 333-334EST
- LMErgoMax™** 333-334XSI
- LMErgoNorm™** 333-334SI



DUAL GRACEYS



LM Dual Gracey™ combines the benefits of a universal curette and a Gracey curette in one. These instruments are designed for the removal of supra and subgingival calculus from all tooth surfaces.

LM™ INNOVATION

Each LM Dual Gracey™ tip has two elliptical cutting edges – a mesial and a distal edge. This design allows for treatment both mesial and distal tooth surfaces. LM Dual Gracays™ are suitable for vertical, horizontal and oblique scaling techniques.



Design:

- Two elliptical cutting edges
- Rounded tip



Effective instrumentation with Dual Gracey

Treatment case by RDH Tatiana Brandt



Syntette™

- A combination of universal curette and Gracey 11/12 and Gracey 13/14 curettes
- For removal of supra and subgingival calculus in premolar and molar areas
- For all tooth surfaces

- LMErgoSense®** 215-216ES
- LMErgoSense®** 215-216EST
- LMErgoMax™** 215-216XSI
- LMErgoNorm™** 215-216SI



- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 215-216SDES
- LMErgoSense®** 215-216SDEST



Mini Syntette™

- Excellent in tight, deep pockets
- Particularly well adapted for use on premolars' and incisors' palatal and lingual surfaces

- LMErgoSense®** 215-216MES
- LMErgoSense®** 215-216MEST
- LMErgoMax™** 215-216MXSI
- LMErgoNorm™** 215-216MSI



- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 215-216MSDES
- LMErgoSense®** 215-216MSDEST



Syntette™ Anterior SD

- A combination of universal curette and Gracey 1/2 and Gracey 7/8 curettes
- For removal of supra and subgingival calculus in premolar and anterior areas
- For all tooth surfaces
- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 275-276SDES
- LMErgoSense®** 275-276SDEST



Mini Syntette™ Anterior SD

- Excellent in narrow, deep pockets
- With sharpen free LM Sharp Diamond™ coated tips

- LMErgoSense®** 275-276MSDES
- LMErgoSense®** 275-276MSDEST

UNIVERSAL CURETTES



Designed for the removal of moderate and heavy calculus from all supra- and subgingival tooth surfaces.


Instrumentation:

For calculus removal place one of the cutting edges of the universal curette against the tooth surface. Tilt the instrument toward the tooth to achieve 70° to 80° angle between the tooth and the blade. Apply lateral pressure and activate the curette by using vertical, diagonal or horizontal pull strokes. Both cutting edges of the same curette can be used by simply tilting the instrument one way or the other.



Design:

- Working end perpendicular to the lower shank
- Two cutting edges (red line)
- Rounded tip

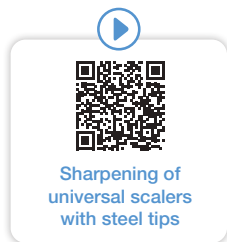

McCALL



McCall 13S-14S

- For incisors and premolars
- Strong blade with rigid shank

LM ErgoSense®	222-223ES
(C)	222-223EST
LM ErgoMax™	222-223XSI
LM ErgoNorm™	222-223SI



Curette LM 15-16

- Alternative for the McCall 13S-14S (222-223)
- For removal of heavy calculus

LM ErgoSense®	225-226ES
(C)	225-226EST
LM ErgoMax™	225-226XSI
LM ErgoNorm™	225-226SI



Mini McCall 13S-14S

- Sharper, slightly rounded tip (not sickle-shaped)
- For deep, narrow pockets

LM ErgoSense®	222-223MES
(C)	222-223MEST
LM ErgoMax™	222-223MXSI
LM ErgoNorm™	222-223MSI



Mini McCall 13S-14S SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	222-223MSDES
(C)	222-223MSDEST



McCall 17-18

- Especially for broad molar surfaces
- Strong, uniform-width working end

LM ErgoSense®	277-278ES
(C)	277-278EST
LM ErgoMax™	277-278XSI
LM ErgoNorm™	277-278SI



McCall 17S-18S

- Excellent in tight interdental spaces
- More delicate blade, working end narrower than in 17-18

LM ErgoSense®	279-280ES
(C)	279-280EST
LM ErgoMax™	279-280XSI
LM ErgoNorm™	279-280SI

LANGER



Langer curettes combine universal curette's working end with Gracey curette's shank angulation. Angulations ease access to difficult areas and both mesial and distal surfaces can be treated with the same instrument end.

Chosen by the Clinician

Effective instrumentation with Langer curettes

Treatment case by RDH Tatiana Brandt



Langer 1/2

- For premolars and molars in the lower jaw

LMErgoSense®	281-282ES
LMErgoSense®	281-282EST
LMErgoMax™	281-282XSI
LMErgoNorm™	281-282SI



Langer 5/6

- For incisors in the upper and lower jaw


LMErgoSense®	285-286ES
LMErgoSense®	285-286EST
LMErgoMax™	285-286XSI
LMErgoNorm™	285-286SI



Langer 3/4

- For premolars and molars in the upper jaw

LMErgoSense®	283-284ES
LMErgoSense®	283-284EST
LMErgoMax™	283-284XSI
LMErgoNorm™	283-284SI

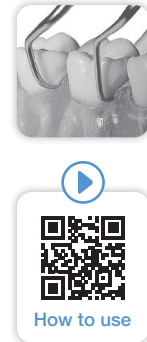

Mini Langers have longer shanks and shorter working ends than standard Langers. Useful in tight interdental spaces. Better access into tight deep pockets and furcations.



Mini Langer 1/2

- For premolars and molars in the lower jaw

LMErgoSense®	281-282MES
LMErgoSense®	281-282MEST
LMErgoMax™	281-282MXSI
LMErgoNorm™	281-282MSI



Mini Langer 3/4

- For premolars and molars in the upper jaw

LMErgoSense®	283-284MES
LMErgoSense®	283-284MEST
LMErgoMax™	283-284MXSI
LMErgoNorm™	283-284MSI



Mini Langer 5/6

- For incisors in the upper and lower jaw

LMErgoSense®	285-286MES
LMErgoSense®	285-286MEST
LMErgoMax™	285-286MXSI
LMErgoNorm™	285-286MSI



GOLDMAN FOX

CE MD
0537



Goldman-Fox 3

- For premolars and molars

LMErgoSense®	233-234ES
☺	233-234EST
LMErgoMax™	233-234XSI
LMErgoNorm™	233-234SI

COLUMBIA

CE MD
0537



Columbia 4L-4R

- For incisors, canines and premolars
- For heavier calculus

LMErgoSense®	219-220ES
☺	219-220EST
LMErgoMax™	219-220XSI
LMErgoNorm™	219-220SI



- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	219-220SDES
☺	219-220SDEST



Columbia 2L-2R

- For incisors and premolars

LMErgoSense®	217-218ES
☺	217-218EST
LMErgoMax™	217-218XSI
LMErgoNorm™	217-218SI



Columbia 13-14

- For molars

LMErgoSense®	237-238ES
☺	237-238EST
LMErgoMax™	237-238XSI
LMErgoNorm™	237-238SI

BARNHART

CE MD
0537



Barnhart 5-6

- For molars
- Shorter lower shank

LMErgoSense®	227-228ES
☺	227-228EST
LMErgoMax™	227-228XSI



- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	227-228SDES
☺	227-228SDEST

DIAMOND FILES

CE MD
0537

Diamond coated files are designed for cleaning deep concave root surfaces in anatomically difficult areas.
Can be used with multidirectional strokes.



Concavity Diamond File SV 5-6

- Diamond coated round curved toe
- For acute furcation lesions

LMErgoSense®	255-256DES
☺	255-256DEST
LMErgoMax™	255-256DXSI
LMErgoNorm™	255-256DSI



Furcation Diamond File SV 7-8

- Diamond coated egg-shaped toe
- For deep furcation lesions

LMErgoSense®	257-258DES
☺	257-258DEST
LMErgoMax™	257-258DXSI
LMErgoNorm™	257-258DSI



Diamond File mes-dist

- Extremely delicate diamond coated instrument for definitive scaling of the root surface
- Easy access with a fine, flexible blade in mesial-distal line angles and developmental grooves


LMErgoSense®	263-264DES
☺	263-264DEST
LMErgoMax™	263-264DXSI



Furcation Diamond File

- Extremely delicate diamond coated instrument for definitive calculus removal in furcations and root depressions
- Easy access with a fine, flexible blade design

LMErgoSense®	261-262DES
☺	261-262DEST
LMErgoMax™	261-262DXSI

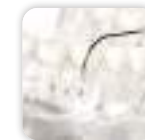

FURCATION INSTRUMENT

CE MD
0537

Well-suited for the cleaning of the furcations of multi-rooted teeth.



Furcator KS

LMErgoSense®	259-260ES
☺	259-260EST
LMErgoMax™	259-260XSI

- Delicate, angled excavator
- Suitable also for the cleaning of narrow approximal concavities and grooves in incisors and premolars
- For cleaning of the furcation roof and inner surfaces, particularly concavations and grooves in them



Handle options:



LMErgoSense® ☺



LMErgoMax™





LMErgoNorm™

HOE SCALERS  
0537



For the removal of supra- and subgingival calculus. Suitable for deep, narrow pockets and concave root surfaces. For both vertical and horizontal techniques.



- Design:**
- Fine hoe-shaped blade with rounded corners
 - Elliptical cutting edge

Sharpening

How to use hoe scalars

Hoe Scaler, anterior

- For the anterior labial and lingual surfaces
- Works well also with lower incisors that are orally inclined

	112-113ES
	112-113EST
	112-113XSI
	112-113SI



Hoe Scaler, lateral


- For the buccal and lingual surfaces of all teeth
- Can also be used for furcations

	134-135ES
	134-135EST
	134-135XSI
	134-135SI



Hoe Scaler, posterior

- For the mesial and distal surfaces of molars
- Can also be used for furcations

	156-157ES
	156-157EST
	156-157XSI
	156-157SI



SURGICAL HANDLE



LM-ErgoHold™ 3

- Excellent grip with elastic silicone coating
- Precise blade control
- Controlled rotation
- Blade position responsiveness maintained
- Controlled movement through tissue
- Fits all blades with fitment No. 3 (blades 10, 10A, 11, 11K, 11P, 12, 12D, 15, 15C and 15S)

	9003
	9003T



I CHOOSE TO WORK SMART.

I love LM Sharp Diamond™ instruments. As a regular user, their superior sharpness, tactile sensitivity, and comfort have truly made a difference in my practice.

With the ergonomic LM-ErgoSense® handle, these instruments feel like natural extensions of my hands, perfectly designed for all forms of scaling. What's truly impressive is their durability. These instruments maintain their original functionality, size, and shape throughout their lifetime, freeing me from frequent maintenance. This saved time is better spent promoting health and crafting beautiful smiles.

The standout feature? Their sharpen free design. Crafted from a durable special metal alloy and enhanced with a protective micro membrane coating, these tools are built to last, enabling me to work efficiently effortlessly. My personal favorite is the LM Dual Gracey™.

It's more than just an instrument; it's my partner in delivering exceptional oral health care.

With LM Sharp Diamond™, I'm not just doing my job - I'm perfecting it. I choose to work smart. I choose to work efficiently. I choose LM Sharp Diamond™.

Dr. Joy Void-Holmes

RDH, BSDH, DHSc
CEO and Founder of Dr. Joy



GRACEY CURETTES



Area specific curettes for the removal of deep subgingival calculus. For the finishing of subgingival root surfaces.



Design:


- Working surface inclined 20° to the lower shank
- One curved cutting edge (red line)
- Rounded cutting tip



Instrumentation:

Place the tip third or half of the cutting edge of a Gracey curette against the tooth surface. The larger outer curve is always the cutting edge. Make sure that the lower shank is parallel with the tooth surface. For calculus removal apply moderate and for root planing light lateral pressure and activate the curette by using vertical, diagonal or horizontal pull strokes.

Area specific indications for Gracey curettes

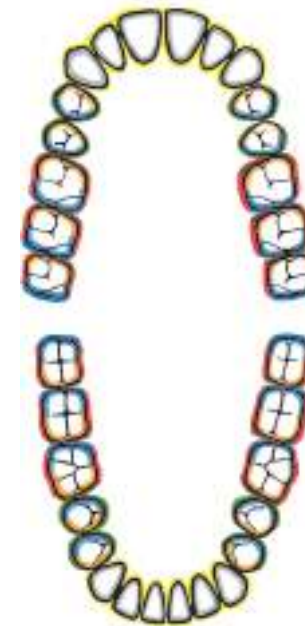

Incisors

○ All surfaces

- Gracey 1/2
- Gracey 5/6
- Gracey 3/4

Premolars and molars

○ Buccal-lingual

- Gracey 7/8
- Gracey 9/10
- Mesial
- Gracey 11/12
- Gracey 15/16

○ Distal

- Gracey 13/14
- Gracey 17/18



Effective instrumentation with Gracey curettes

- Mini Gracey 1/2
- Mini Gracey 11/12
- Macro Gracey 13/14
- Mini Sickle

Treatment case by RDH Tatiana Brandt



Handle options:



LMErgoSense®



LMErgoSense®



LMErgoMax™



LMErgoNorm™

STANDARD AND MINI GRACEY



Mini Gracey has longer lower shank and smaller working end than in the standard model. Functions well in deep narrow pockets.



Gracey 1/2

- For incisors and premolars

LMErgoSense®	201-202ES
☞	201-202EST
LMErgoMax™	201-202XSI
LMErgoNorm™	201-202SI



- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	201-202SDES
☞	201-202SDEST



Mini Gracey 1/2

- For incisors and premolars

LMErgoSense®	201-202MES
☞	201-202MEST
LMErgoMax™	201-202MXSI
LMErgoNorm™	201-202MSI



- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	201-202MSDES
☞	201-202MSDEST



Gracey 3/4

- For incisors and premolars

LMErgoSense®	203-204ES
☞	203-204EST
LMErgoMax™	203-204XSI
LMErgoNorm™	203-204SI



Gracey P3-P4

- Short shank makes the instrument more rigid

LMErgoSense®	266-267ES
☞	266-267EST
LMErgoMax™	266-267XSI
LMErgoNorm™	266-267SI



Gracey 5/6

- For incisors and premolars
- Especially useful in proximal surfaces

LM ErgoSense®	205-206ES
LM ErgoSense®	205-206EST
LM ErgoMax™	205-206XSI
LM ErgoNorm™	205-206SI



Gracey 7/8

- For premolars and molars
- For buccal, mesial and lingual surfaces

LM ErgoSense®	207-208ES
LM ErgoSense®	207-208EST
LM ErgoMax™	207-208XSI
LM ErgoNorm™	207-208SI



Gracey 7/8 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	207-208SDES
LM ErgoSense®	207-208SDEST



Mini Gracey 7/8

- For premolars and molars
- For buccal, mesial and lingual surfaces

LM ErgoSense®	207-208MES
LM ErgoSense®	207-208MEST
LM ErgoMax™	207-208MXSI
LM ErgoNorm™	207-208MSI



Mini Gracey 7/8 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	207-208MSDES
LM ErgoSense®	207-208MSDEST



Gracey 9/10

- For molars' buccal and lingual surfaces using horizontal technique

LM ErgoSense®	209-210ES
LM ErgoSense®	209-210EST
LM ErgoMax™	209-210XSI
LM ErgoNorm™	209-210SI



Gracey 9/10 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	209-210SDES
LM ErgoSense®	209-210SDEST



Gracey 11/12

- For premolars' and molars' mesial surfaces

LM ErgoSense®	211-212ES
(C)	211-212EST
LM ErgoMax™	211-212XSI
LM ErgoNorm™	211-212SI



Gracey 11/12 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	211-212SDES
(C)	211-212SDEST



Mini Gracey 11/12

- For premolars' and molars' mesial surfaces

LM ErgoSense®	211-212MES
(C)	211-212MEST
LM ErgoMax™	211-212MXSI
LM ErgoNorm™	211-212MSI



Mini Gracey 11/12 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	211-212MSDES
(C)	211-212MSDEST



Gracey 13/14

- For premolars' and molars' distal surfaces

LM ErgoSense®	213-214ES
(C)	213-214EST
LM ErgoMax™	213-214XSI
LM ErgoNorm™	213-214SI



Gracey 13/14 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	213-214SDES
(C)	213-214SDEST



Mini Gracey 13/14

- For premolars' and molars' distal surfaces

LM ErgoSense®	213-214MES
(C)	213-214MEST
LM ErgoMax™	213-214MXSI
LM ErgoNorm™	213-214MSI



Mini Gracey 13/14 SD

- With sharpen free LM Sharp Diamond™ coated tips

LM ErgoSense®	213-214MSDES
(C)	213-214MSDEST



Gracey 15/16

- For molars' mesial surfaces
- Blade shank more sharply angulated than in Gracey 11/12 for better access

LMErgoSense®	245-246ES
◀	245-246EST
LMErgoMax™	245-246XSI
LMErgoNorm™	245-246SI



Mini Gracey 15/16

- Blade shank more sharply angulated than in Mini Gracey 11/12 for better access

LMErgoSense®	245-246MES
◀	245-246MEST
LMErgoMax™	245-246MXSI
LMErgoNorm™	245-246MSI



Gracey 15/16 SD

- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	245-246SDES
◀	245-246SDEST



Mini Gracey 15/16 SD

- With sharpen free LMSharp Diamond™ coated tips


LMErgoSense®	245-246MSDES
◀	245-246MSDEST



Gracey 17/18

- For molars' distal surfaces
- Blade shank more sharply angulated than in Gracey 13/14 for better access

LMErgoSense®	247-248ES
◀	247-248EST
LMErgoMax™	247-248XSI
LMErgoNorm™	247-248SI

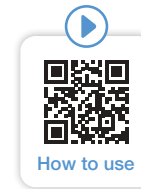

Mini Gracey 17/18

- Blade shank more sharply angulated than in Mini Gracey 13/14 for better access

LMErgoSense®	247-248MES
◀	247-248MEST
LMErgoMax™	247-248MXSI
LMErgoNorm™	247-248MSI



Gracey 17/18 SD

- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	247-248SDES
◀	247-248SDEST



Mini Gracey 17/18 SD

- With sharpen free LM Sharp Diamond™ coated tips

LMErgoSense®	247-248MSDES
◀	247-248MSDEST





MACRO GRACEY



Standard model's long working end combined with mini model's longer lower shank. Excellent reach interdentally in deep pockets.



Macro Gracey 5/6

- For incisors and premolars
- Especially useful in proximal surfaces

 205-206AFES	
 205-206AFEST	
 205-206AFXSI	
 205-206AFSI	





Macro Gracey 7/8

- For premolars and molars
- For buccal, mesial and lingual surfaces

 207-208AFES	
 207-208AFEST	
 207-208AFXSI	
 207-208AFSI	





Macro Gracey 11/12

- For premolars' and molars' mesial surfaces

 211-212AFES	
 211-212AFEST	
 211-212AFXSI	
 211-212AFSI	

Macro Gracey 13/14

- For premolars' and molars' distal surfaces

 213-214AFES	
 213-214AFEST	
 213-214AFXSI	
 213-214AFSI	

RIGID GRACEY



Blade shank thicker and more rigid, and working end slightly wider than in the standard model. Permits the use of greater pressure.

Rigid Gracey 1/2

- For incisors and premolars

 201-202RES	
 201-202REST	
 201-202RXSI	
 201-202RSI	

Rigid Gracey 7/8

- For premolars and molars
- For buccal, mesial and lingual surfaces

 207-208RES	
 207-208REST	
 207-208RXSI	
 207-208RSI	

Rigid Gracey 11/12

- For premolars' and molars' mesial surfaces

 211-212RES	
 211-212REST	
 211-212RXSI	
 211-212RSI	


Rigid Gracey 13/14

- For premolars' and molars' distal surfaces


 213-214RES	
 213-214REST	
 213-214RXSI	
 213-214RSI	

Rigid Gracey 15/16

- For molars' mesial surfaces
- Blade shank more sharply angulated than in Gracey 11/12 for better access

 245-246RES	
 245-246REST	
 245-246RXSI	
 245-246RSI	

The types of Gracey curettes:

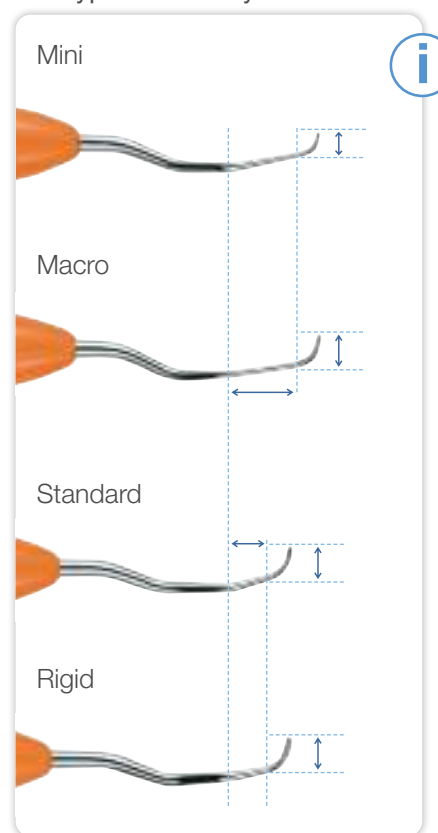

MESIAL-DISTAL GRACEY



Combining a mesial and a distal working end in one instrument. Useful when treating one quadrant at a time.





All mesial-distal Gracays have code rings that make it easy to identify which working end is mesial and which distal.

Distal





Mesial





Gracey 11/14

	211-214ES
	211-214EST
	211-214XSI
	211-214SI





Gracey 12/13

	212-213ES
	212-213EST
	212-213XSI
	212-213SI

Gracey 15/18

	245-248ES
	245-248EST
	245-248XSI
	245-248SI





Gracey 16/17

	246-247ES
	246-247EST
	246-247XSI
	246-247SI





MINI MESIAL-DISTAL GRACEY

Mini models have longer lower shank and smaller working end than the standard model. Functions well in deep narrow pockets.





Mini Gracey 11/14

	211-214MES
	211-214MEST
	211-214MXSI
	211-214MSI





Mini Gracey 12/13

	212-213MES
	212-213MEST
	212-213MXSI
	212-213MSI

Mini Gracey 15/18

	245-248MES
	245-248MEST
	245-248MXSI
	245-248MSI





Mini Gracey 16/17

	246-247MES
	246-247MEST
	246-247MXSI
	246-247MSI





MACRO MESIAL-DISTAL GRACEY

Macro models have standard model's long working end combined with mini model's longer lower shank. Excellent reach interdentally in deep pockets.

Macro Gracey 11/14


	211-214AFES
	211-214AFEST
	211-214AFXSI
	211-214AFSI

Macro Gracey 12/13

	212-213AFES
	212-213AFEST
	212-213AFXSI
	212-213AFSI



Handle options:

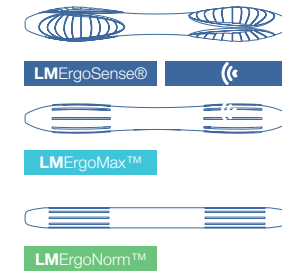

EXCHANGEABLE TITANIUM TIPS FOR IMPLANTS



LMErgoMix™


A great instrument series for implant maintenance. The LMErgoMix™ structure ensures excellent tactile sensitivity and optimal rigidity.

LM-ErgoMix™ hand instrument series offers the dental team a specially economical, ecological and ergonomic alternative. Whenever there is a need for a tip replacement the lock grip is opened, a new tip is inserted and the lock grip is closed again. No extra tools are needed.

Use LM-ErgoMix™ instruments in the normal way, clean and sterilize it like any other instrument. The sealed construction guarantees high level hygiene.

Benefits of retipping:

- Always sharp instruments
- Minimizes waste
- Enables personalized tip combinations.

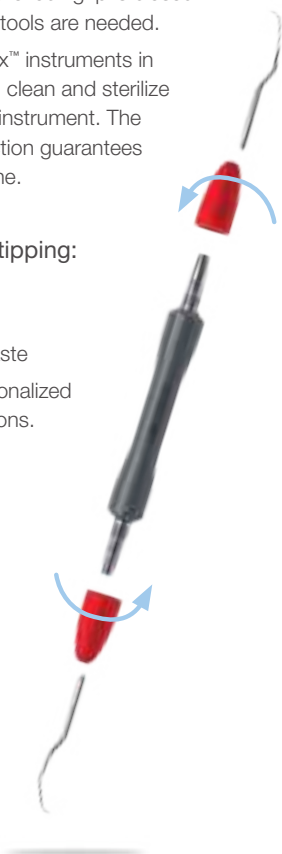

Titanium tips for sensitive but still effective implant maintenance.



Implant Mini Universal curette

- Delicate working end
- Especially for tight interdental spaces

LMErgoMix™ 283-284MTIEM



Implant Mini Gracey 1/2

- For incisors and premolars

LMErgoMix™ 201-202MTIEM



Implant Mini Gracey 11/12

- For premolars' and molars' mesial surfaces

LMErgoMix™ 211-212MTIEM



Implant Mini Gracey 13/14

- For premolars' and molars' distal surfaces

LMErgoMix™ 213-214MTIEM

LM-ErgoMix™ spare parts

Spare handle	
LM-ErgoMix™ handle dark grey	9011EM
Spare lock grips	
Orange	9020EM
Light blue	9021EM
Grey	9022EM
Yellow	9023EM
Dark grey	9024EM
Red	9025EM
Spare titanium tips	
Implant Mini Universal Curette L	283MTIEM
Implant Mini Universal Curette R	284MTIEM
Implant Mini Gracey 1	201MTIEM
Implant Mini Gracey 2	202MTIEM
Implant Mini Gracey 11	211MTIEM
Implant Mini Gracey 12	212MTIEM
Implant Mini Gracey 13	213MTIEM
Implant Mini Gracey 14	214MTIEM

PERIODONTAL INSTRUMENT SETS

CE MD
0537



LM-PerioSet 1

- LM**ErgoSense® 6652ES
- LM**ErgoMax™ 6652XSI

- Includes:**
- Mini Sickle
 - Curette McCall 17-18
 - Gracey 1/2
 - Gracey 11/12
 - Gracey 13/14
 - LM-Servo™ 5 cassette

LM-PerioSet 2

- LM**ErgoSense® 6682ES
- LM**ErgoMax™ 6682XSI

- Includes:**
- Mini Sickle
 - Curette McCall 17-18
 - Gracey 1/2
 - Gracey 9/10
 - Gracey 11/12
 - Gracey 13/14
 - Mini Gracey 15/16
 - Mini Gracey 17/18
 - LM-Servo™ 8 cassette

LM-PerioKit

- LM**ErgoSense® 6842ES

- Includes:**
- Implant Mini Universal curette
 - Concavity File SV 5-6
 - Mini Gracey 13-14
 - Mini Gracey 11-12
 - Mini Gracey 7-8
 - McCall 17S-18S
 - North Carolina Probe
 - Furcation Probe
 - LM-Servo™ 8 cassette

”

The LM-PerioKit by StyleItaliano™ contains the essentials for the best possible therapy.

Instruments have been selected in the spirit of StyleItaliano™ whose philosophy is to make treatment feasible, teachable, repeatable. Periodontal therapy too can follow this approach. It's my feeling that, in the best interest of his patients, every practitioner should be able to properly execute non surgical therapy, crown lengthening procedures, free gingival grafts, simpler surgical root coverages.”

Dr. Giano Ricci


MD, DDS, MScD



SPOON EXCAVATORS



For the removal of softened carious dentin and temporary fillings.

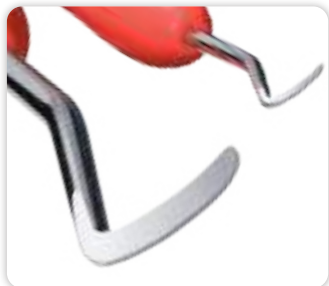

Design:

- Spoon-shaped working end
- Rounded tip
- Long cutting edge, which can easily treat broader areas



Handle options:



Spoon Excavator, 1.0 mm

- LM ErgoSense® 522-523ES
- CE 522-523EST
- LM ErgoMax™ 522-523XSI
- LM ErgoNorm™ 522-523SI



Spoon Excavator, 1.5 mm

- LM ErgoSense® 533-534ES
- CE 533-534EST
- LM ErgoMax™ 533-534XSI
- LM ErgoNorm™ 533-534SI

ROUND EXCAVATORS



Cariology:

- For the removal of softened, carious dentin
- For the removal of temporary fillings

Surgery:

- To sever the periodontal ligament prior to extraction
- For the removal of granulation tissue after extraction




Design:

- Round working end
- Curved cutting edge



Round Excavator, 1.0 mm

- LM ErgoSense® 61-62ES
- CE 61-62EST
- LM ErgoMax™ 61-62XSI
- LM ErgoNorm™ 61-62SI

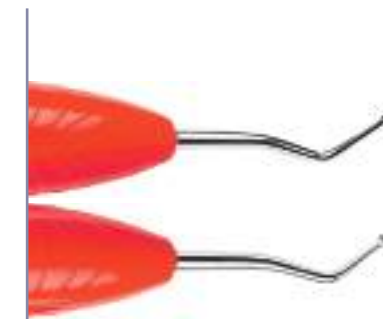

Round Excavator, 1.5 mm

- LM ErgoSense® 63-64ES
- CE 63-64EST
- LM ErgoMax™ 63-64XSI
- LM ErgoNorm™ 63-64SI



Round Excavator, 2.0 mm

- LM ErgoSense® 65-66ES
- CE 65-66EST
- LM ErgoMax™ 65-66XSI
- LM ErgoNorm™ 65-66SI

TRIPLE-ANGLED BACK ACTION EXCAVATORS



Cariology:

- For removing carious dentin
- For undermining enamel during tissue-saving preparations
- For the horizontal finishing of outer cavity margins

Periodontics:

- For the removal of calculus from hard to reach areas such as the lower anterior lingual surfaces
- For the depuration of concave root grooves and furcations

Surgery:

- For the removal of granulation tissue in connection with apical resections, flap operations, or extractions



Design:

- Round working end and curved cutting edge
- Shank angulations ease treatments on hard to reach buccal, distal, and lingual surfaces



The blade shank is triple-angled, allowing for "back action" functions



Triple-angled Back Action Excavator, 1.0 mm

	611-621ES
	611-621EST
	611-621XSI
	611-621SI



Triple-angled Back Action Excavator, 1.5 mm

	631-641ES
	631-641EST
	631-641XSI
	631-641SI



Triple-angled Back Action Excavator, 2.0 mm

	651-661ES
	651-661EST
	651-661XSI
	651-661SI



Triple-angled Back Action Excavator, 3.5 mm

- For the horizontal finishing of outer cavity margins

	671-681ES
	671-681EST
	671-681XSI
	671-681SI

GINGIVAL MARGIN TRIMMERS



For finishing enamel margins in the proximal area as well as for making plastic, amalgam and ceramic fillings.



Gingival Margin Trimmer, 1.2 mm, distal

	121-122ES
	121-122EST
	121-122XSI
	121-122SI



Gingival Margin Trimmer, 1.2 mm, mesial

	125-126ES
	125-126EST
	125-126XSI
	125-126SI

SPATULA



Spatula

- LMErgoSense®** 450-460ES
- ☺** 450-460EST
- LMErgoMax™** 450-460XSI
- LMErgoNorm™** 450-460SI

CORD PACKERS



For packing retraction cord into the gingival pocket.



The bent blade enables getting around the tooth as well as the packing of the cord using vertical motions.



Cord Packer, 2.0 mm

- LMErgoSense®** 490-491ES
- ☺** 490-491EST
- LMErgoMax™** 490-491XSI
- LMErgoNorm™** 490-491SI



Cord Packer, 2.6 mm

- LMErgoSense®** 492-493ES
- ☺** 492-493EST
- LMErgoMax™** 492-493XSI
- LMErgoNorm™** 492-493SI

BASE CARRIERS



For carrying base or thin liner to dentin surfaces.



Base Carrier, rigid

- LMErgoSense®** 29H-30HES
- ☺** 29H-30HEST
- LMErgoMax™** 29H-30HXS
- LMErgoNorm™** 29H-30HSI



Base Carrier – Condenser, malleable tip

- Tip (30) can be bent into the desired form, allowing for easier access to hard to reach areas

- LMErgoSense®** 30-34ES
- ☺** 30-34EST
- LMErgoMax™** 30-34XSI
- LMErgoNorm™** 30-34SI



Base Carrier - Condenser, rigid

- LMErgoSense®** 30H-34ES
- ☺** 30H-34EST
- LMErgoMax™** 30H-34XSI
- LMErgoNorm™** 30H-34SI

PLUGGERS

For packing all types of filling materials.



Plugger, cylindrical,
1.2-2.0 mm

- LMErgoSense®** 350-380ES
- LMErgoSense®** 350-380EST
- LMErgoMax™** 350-380XSI
- LMErgoNorm™** 350-380SI



Plugger, cylindrical,
1.8-2.5 mm

- LMErgoSense®** 370-390ES
- LMErgoSense®** 370-390EST
- LMErgoMax™** 370-390XSI
- LMErgoNorm™** 370-390SI

Plugger, cylindrical,
1.0-1.5 mm

- LMErgoSense®** 340-360ES
- LMErgoSense®** 340-360EST
- LMErgoMax™** 340-360XSI
- LMErgoNorm™** 340-360SI

BURNISHERS

For burnishing, modelling and finishing of fillings.



Burnisher

- For burnishing fillings at cavity margins
- For packing and shaping of composite fillings
- For packing of inlay-/onlay fillings

- LMErgoSense®** 31-32ES
- LMErgoSense®** 31-32EST
- LMErgoMax™** 31-32XSI
- LMErgoNorm™** 31-32SI



Burnisher DD

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 31-32DDES
- LMErgoSense®** 31-32DDEST



Burnisher – Plugger,
ball-ended, 1.5-2.5 mm

- For plugging filling into cavity
- For polishing
- For plugging and shaping of composite fillings

- LMErgoSense®** 471-473ES
- LMErgoSense®** 471-473EST
- LMErgoMax™** 471-473XSI
- LMErgoNorm™** 471-473SI



Burnisher – Plugger DD,
ball-ended, 1.5-2.5 mm

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 471-473DDES
- LMErgoSense®** 471-473DDEST



Burnisher 27–29

- For shaping of occlusal surfaces
- For burnishing of fillings at cavity margins

- LMErgoSense®** 476-477ES
- LMErgoSense®** 476-477EST
- LMErgoMax™** 476-477XSI
- LMErgoNorm™** 476-477SI



Burnisher 21

- For shaping and packing of composite fillings
- For shaping occlusal surfaces before hardening
- Pits, grooves, lateral ridges

- LMErgoSense®** 476-477ES
- LMErgoSense®** 476-477EST
- LMErgoMax™** 476-477XSI
- LMErgoNorm™** 476-477SI

LM-ARTE™ FOR ESTHETIC RESTORATIONS



LM-Arte™ is a set of innovative instruments for esthetic restorations and it is designed especially for composite layering.



CO-OPERATION WITH STYLEITALIANO™

LM-Arte™ instruments have been developed in collaboration with Styleitaliano™, a group of passionate clinicians and researchers specialized in esthetic restorative dentistry. The unique concept of LM-Arte™ instruments has been a breakthrough in esthetic dentistry.

Read more on pages 10-11.



REPLICATE EASILY

I can make modeling within seconds with simple steps by using LM-Arte™ Replica Posterior instrument and technique.

Dr. Osama Shaalan
Gold Member of Styleitaliano™



LM-Arte™ Replica Posterior

- For easy occlusal modeling of posterior teeth in direct composite restorations
- Occlusal anatomy of single tooth is guided by buccal cusps using replica technique: adapt, shape, replicate
- Concave and convex tip profiles

LMErgoSense® 448-449ES

LMErgoMax™ 448-449XSI

LM-Arte™ Solo instruments are ideal for simplifying your daily cases: with just two instruments, you can perform most of your high-end anterior and posterior restorations.



LM-Arte™ Solo Anterior

- For high-end esthetic anterior restorations
- Rounded and straight spatulas are designed for modeling anterior teeth, especially labial surfaces

LMErgoSense® 431-433ES

LMErgoSense® 431-433EST

LMErgoMax™ 431-433XSI


LM-Arte™ Solo Posterior

- For high-end posterior composite restorations
- The plugger end is designed for bulk filling and composite material modeling
- The sharp end enables contouring of the occlusal surface anatomy

LMErgoSense® 434-438ES

LMErgoSense® 434-438EST

LMErgoMax™ 434-438XSI

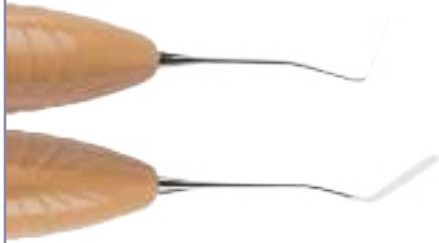


LM-Arte™ Applica

- A very thin and flexible spatula for transporting and modeling of the composite



- LMErgoSense®** 46-49ES
- LMErgoSense®** 46-49EST
- LMErgoMax™** 46-49XSI

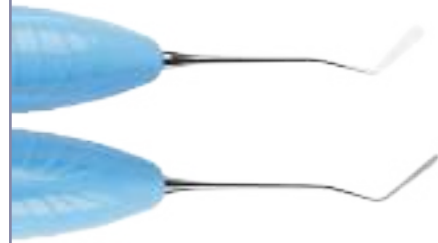

LM-Arte™ Modella

- A very thin and flexible spatula for transporting and modeling of the composite



- LMErgoSense®** 442-443ES
- LMErgoSense®** 442-443EST
- LMErgoMax™** 442-443XSI



LM-Arte™ Condensa

- A long and round plunger especially for modeling composite



- LMErgoSense®** 488-489ES
- LMErgoSense®** 488-489EST
- LMErgoMax™** 488-489XSI



LM-Arte™ Fissura

- Sharp pointed tips for modeling of occlusal surface structures



- LMErgoSense®** 481-487ES
- LMErgoSense®** 481-487EST
- LMErgoMax™** 481-487XSI



NEW!



LM-Arte™ Applica DD

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 46-49DDES
- LMErgoSense®** 46-49DDEST

NEW!



LM-Arte™ Modella DD

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 442-443DDES
- LMErgoSense®** 442-443DDEST

NEW!



LM-Arte™ Condensa DD

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 488-489DDES
- LMErgoSense®** 488-489DDEST

NEW!



LM-Arte™ Fissura DD

- With LM Dark Diamond™ non-stick coated tips

- LMErgoSense®** 481-487DDES
- LMErgoSense®** 481-487DDEST



LM-Arte™ Eccesso

- For the removal of composite and bonding residue

LMErgoSense®	307-308ES
◀	307-308EST
LMErgoMax™	307-308XSI



LM-Arte™ Misura

- For horizontal and vertical measuring of thickness of composite layers on the restorations

LMErgoSense®	496-497ES
◀	496-497EST
LMErgoMax™	496-497XSI



LM-Arte™ Cusp Misura

- Unique instrument for estimating cusp thickness and height for evaluating the need for cusp preparation when performing adhesive restorations.

LMErgoSense®	468-469ES
◀	468-469EST
LMErgoMax™	468-469XSI



LM-Arte™ Applica Twist

- Thin and narrow spatula with bended angles for modeling difficult to reach areas such as approximal walls or the marginal ridge in matrix

LMErgoSense®	464-494ES
◀	464-494EST
LMErgoMax™	464-494XSI



LM-Arte™ Posterior Misura

- For horizontal and vertical measuring of thickness of composite layers on the posterior restorations
- The fork shaped tip allows measuring and modeling of the composite wall with the presence of a matrix

LMErgoSense®	498-499ES
◀	498-499EST
LMErgoMax™	498-499XSI



LM-Arte™ Eccesso DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	307-308DDES
◀	307-308DDEST

NEW!




LM-Gengiva™ gingival retracting kit

Includes:

- LM-Gengiva™ / LM-MultiHolder™ PK III hand instrument (see also page 53)
- LM-Gingival Retractors (blue) for establishing optimum working conditions for procedures at the gum line. Assortment package of three sizes include circa 90 pieces, 30 pcs/size): S for incisors, M for premolars and XL for molars
- LM-MultiLever™ for the removal of disposable gingival retractors

LMErgoSense®	7556KES
LMErgoMax™	7556KXSI

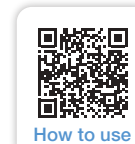

LM-ARTE™ INSTRUMENT SETS



The non-stick LM Dark Diamond™ coated LM-Arte™ instruments feature exclusive coating for smooth and precise non-stick placement and sculpting of composite resins.

What do the StyleItaliano™ experts think? Watch the video!



LM-Arte™ Set

LMErgoSense® 6840ES

LMErgoMax™ 6840XSI

Includes:

- LM-Arte™ Eccesso
- LM-Arte™ Misura
- LM-Arte™ Fissura
- LM-Arte™ Applica
- LM-Arte™ Condensa
- LM-Servo™ 5 cassette

LM-Arte™ Posterior Set

LMErgoSense® 6844ES

LMErgoMax™ 6844XSI

Includes:

- LM-Arte™ Posterior Misura
- LM-Arte™ Cusp Misura
- LM-Arte™ Fissura
- LM-Arte™ Modella
- LM-Arte™ Condensa
- LM-Servo™ 5 cassette

LM-Arte™ Complete Set

LMErgoSense® 6848ES

LMErgoMax™ 6848XSI

Includes:

- LM-Arte™ Eccesso
- LM-Arte™ Misura
- LM-Arte™ Applica
- LM-Arte™ Posterior Misura
- LM-Arte™ Cusp Misura
- LM-Arte™ Fissura
- LM-Arte™ Modella
- LM-Arte™ Condensa
- LM-Servo™ 8 cassette

NEW!

LM-Arte™ Dark Diamond set

LMErgoSense® 6845ES

- Instruments with LM Dark Diamond™ non-stick coated tips

Includes:

- LM-Arte™ Eccesso DD
- LM-Arte™ Modella DD
- LM-Arte™ Fissura DD
- LM-Arte™ Applica DD
- LM-Arte™ Condensa DD
- LM-Servo™ 5 cassette

FOR COMPOSITES **MD** **CE**



Modeller MH

- For packing large, approximal fillings in several layers
- Cylinder shaped working ends won't get stuck into the composite while light curing
- For pressing the matrix band in place while curing

LMErgoSense®	708-709ES
☺	708-709EST
LMErgoMax™	708-709XSI



Modeller MH DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	708-709DDES
☺	708-709DDEST



Modeller MH Mini

- For packing small or medium approximal fillings in several layers
- Cylinder shaped working ends won't get stuck into the composite while light curing
- For pressing the matrix band in place while curing

LMErgoSense®	707-708ES
☺	707-708EST
LMErgoMax™	707-708XSI



Modeller MH Mini DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	707-708DDES
☺	707-708DDEST



Flow Instrument IO

- For easy handling and shaping of flowable composites

LMErgoSense®	484-485ES
☺	484-485EST
LMErgoMax™	484-485XSI



Form Condenser, 1.5-2.5 mm

- For shaping and packing of large filling surfaces

LMErgoSense®	702-703ES
☺	702-703EST
LMErgoMax™	702-703XSI
LMErgoNorm™	702-703SI



Composite Spatula 6

- For placing and contouring composite

LMErgoSense®	446-447ES
☺	446-447EST
LMErgoMax™	446-447XSI
LMErgoNorm™	446-447SI



Handle options:





Packer-Modeller IO p

- Combination instrument for premolars


LMErgoSense®	482-702ES
☞	482-702EST
LMErgoMax™	482-702XSI
LMErgoNorm™	482-702SI



Packer-Modeller IO m

- Combination instrument for molars

LMErgoSense®	483-703ES
☞	483-703EST
LMErgoMax™	483-703XSI
LMErgoNorm™	483-703SI

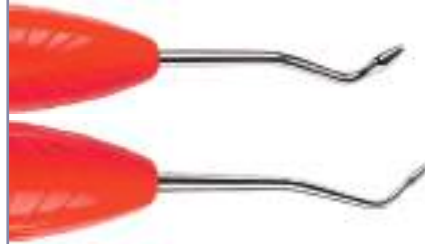

Packer-Modeller IO

- Conical working ends, which are exceptionally well-suited e.g. to shaping of occlusal surfaces and to packing of all kinds of plastics
- Composites do not stick to the highly polished blades



LMErgoSense®	482-483ES
☞	482-483EST
LMErgoMax™	482-483XSI
LMErgoNorm™	482-483SI



Packer-Modeller IO DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	482-483DDES
☞	482-483DDEST



Plastic Instrument, anterior

- For placing, packing and contouring composite

LMErgoSense®	43-47ES
☞	43-47EST
LMErgoMax™	43-47XSI
LMErgoNorm™	43-47SI



Plastic instrument, anterior DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	43-47DDES
☞	43-47DDEST



Plastic Instrument, posterior

- For placing, packing and contouring composite

LMErgoSense®	48-702ES
☞	48-702EST
LMErgoMax™	48-702XSI
LMErgoNorm™	48-702SI



Plastic instrument, posterior DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	48-702DDES
☞	48-702DDEST

CARVERS



For shaping of filling material to fit the anatomy of tooth, and for removal of soft excess.



Mini Spatula

- For placing and contouring composite

LMErgoSense®	444-445ES
	444-445EST
LMErgoMax™	444-445XSI
LMErgoNorm™	444-445SI



Mini Spatula DD

- With LM Dark Diamond™ non-stick coated tips

LMErgoSense®	444-445DDES
	444-445DDEST



1/2 Hollenback

- For proximal and cervical surfaces

LMErgoSense®	741-751ES
	741-751EST
LMErgoMax™	741-751XSI
LMErgoNorm™	741-751SI



Carver Maehlum

- For removal of soft excess from gingival margins (71)
- For anatomical shaping of occlusal surfaces (73)

LMErgoSense®	71-73ES
	71-73EST
LMErgoMax™	71-73XSI
LMErgoNorm™	71-73SI



Carver Nyström I

- For proximal and occlusal surfaces

LMErgoSense®	74-75ES
	74-75EST
LMErgoMax™	74-75XSI
LMErgoNorm™	74-75SI



Carver Nyström III

- For interdental spaces, contact surfaces

LMErgoSense®	76-77ES
	76-77EST
LMErgoMax™	76-77XSI
LMErgoNorm™	76-77SI



Universal Carver, occlusal-proximal

- For removal of soft excess from gingival margins (71)
- For anatomical shaping of occlusal surfaces (72)

LMErgoSense®	71-72ES
	71-72EST
LMErgoMax™	71-72XSI
LMErgoNorm™	71-72SI



Cleoid-Discoïd

- For anatomical shaping of occlusal surfaces of small fillings

LMErgoSense®	731-732ES
	731-732EST
LMErgoMax™	731-732XSI
LMErgoNorm™	731-732SI

EXCESS SCALERS

MD CE



Excess Scaler

- For the removal of adhesive, cement, and composite excesses

LMErgoSense® 785-795ES

 785-795EST

LMErgoMax™ 785-795XSI

LMErgoNorm™ 785-795SI



Mini Excess Scaler

- An effective excess remover in extremely tight interdental spaces

LMErgoSense® 784-794ES

 784-794EST

LMErgoMax™ 784-794XSI

LMErgoNorm™ 784-794SI

FILLING REMOVER

MD CE



Filling Remover

- For the removal of temporary fillings and crowns
- For bridge removal during fitting

LMErgoSense® 764-766ES

 764-766EST

LMErgoMax™ 764-766XSI

LMErgoNorm™ 764-766SI

CROWN OFF

MD CE



A groove is first burred from the lower edge of the crown up to the occlusal surface. The groove should be so deep that it reaches the cement. The Crown Off working end is placed longitudinally into the groove and turned as if spreading the groove. The blade is turned until the crown loosens and lifts up. The blade stays firmly in the groove and it is easy to control even when lot of force is used.



Crown Off

- Blade double angled with a 45° twist
- Crown removal instrument

LMErgoSense® 767-768ES

 767-768EST

LMErgoMax™ 767-768XSI

INSTRUMENT SET

MD CE



LM-CompoSet

LMErgoSense® 6657ES

LMErgoMax™ 6657XSI

Includes:


- Mini Spatula
- Plastic instrument, anterior
- Form Condenser
- Packer-Modeller IO
- Mini Excess Scaler
- LM-Servo™ 5 cassette

LM-MULTIHOLDER™ SYSTEM **MD** **CE**

LM-MultiHolder™ instruments and different kinds of replaceable plastic working ends – a real multi-function instrument!

Angle options allow for easy access to even the more complex areas. Optimal combination of steel and plastic: the steel shank is rigid and only the working end will flex, when required.

LM™ INNOVATION

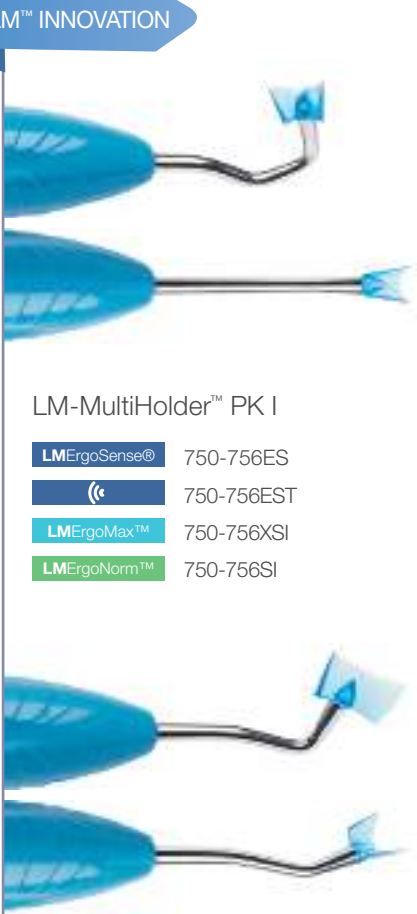

LM-MultiHolder™ PK I

- LM**ErgoSense® 750-756ES
- (C)** 750-756EST
- LM**ErgoMax™ 750-756XSI
- LM**ErgoNorm™ 750-756SI



LM-MultiHolder™ PK II

- LM**ErgoSense® 752-754ES
- (C)** 752-754EST
- LM**ErgoMax™ 752-754XSI
- LM**ErgoNorm™ 752-754SI

LM-MultiHolder™ PK III

• See page 47 for LM-Gengiva™

- LM**ErgoSense® 752-756ES
- (C)** 752-756EST
- LM**ErgoMax™ 752-756XSI

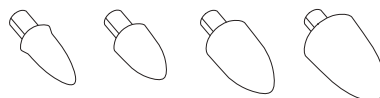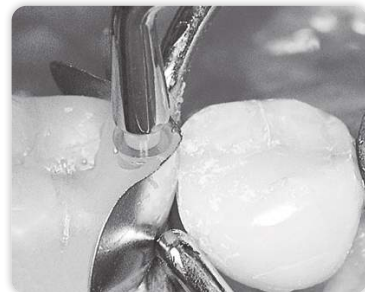


LM-Contact Formers

Transparent tip for making correctly formed and tight approximal contact.

- Suitable for use with all matrices
- Will not stick to composite fillings
- Disposable

Four sizes corresponding to the size of the cavity, circa 60 pcs/package:
 7622 XS
 7624 S
 7626 M
 7628 L
 7620 Assortment package circa 120 pcs (30 pcs/size)

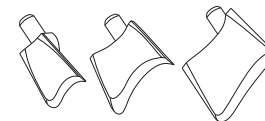


LM-Gingival Retractors (blue)

For establishing optimum working conditions for procedures at the gum line.

- Disposable

Three sizes, circa 60 pcs/package:
 7502 S for incisors
 7504 M for premolars
 7506 XL for molars
 7500 Assortment package circa 90 pcs (30 pcs/size)

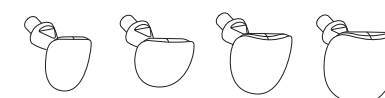

LM-Cervical Matrices

For making anatomically correct fillings at the gum line.

- Protects from moisture
- Prevents air bubbles
- Disposable

Four models for different size and shape, circa 60 pcs/package:
 7602 for incisors
 7604 for canines and premolars
 7606 for small molars
 7608 for large molars
 7600 Assortment package circa 120 pcs (30 pcs/size)

LM-MultiHolder™ IntroKit

Introkit contains LM-MultiHolder™ instruments PK I and PK II, LM-MultiLever™ and working end sample package.
 (Incl 44 working ends; 4 for each size of each working end type)

- 7555KES
- 7555KXSI
- Separate working end sample package 7555



LM-MultiLever™ 7550

- For the removal of plastic working end

ROOT CANAL EXPLORER



The instrument's long, sharp tips facilitate the location and opening of root canals to provide access for proper preparation instruments.

Important! Only use the instruments for the purpose which they are intended!



Root Canal Explorer

- LMErgoSense®** 27-28ES
- LMErgoSense®** 27-28EST
- LMErgoMax™** 27-28XSI
- LMErgoNorm™** 27-28SI

ROOT CANAL EXCAVATORS



Root canal excavators have longer lower shanks and sharper blade angles than round excavators.

- Extended and correctly angled shanks as well as highly effective blade angles allow for the preparation of all types of pulp cavities
- Effective for the cleaning of pulp chambers and location of root canals as well as the removal of temporary fillings



Handle options:



LMErgoSense®



LMErgoMax™



LMErgoNorm™



Root Canal Excavator, 1.0 mm

- LMErgoSense®** 612-622ES
- LMErgoSense®** 612-622EST
- LMErgoMax™** 612-622XSI
- LMErgoNorm™** 612-622SI



Root Canal Excavator, 1.5 mm

- LMErgoSense®** 632-642ES
- LMErgoSense®** 632-642EST
- LMErgoMax™** 632-642XSI
- LMErgoNorm™** 632-642SI



Root Canal Excavator, 2.0 mm

- LMErgoSense®** 652-662ES
- LMErgoSense®** 652-662EST
- LMErgoMax™** 652-662XSI
- LMErgoNorm™** 652-662SI

ENDO SPREADERS




For the compacting of root canal filling materials during lateral condensation.



Endo Spreader 0.4 mm

- LMErgoSense®** 24ES
- LMErgoSense®** 24EST
- LMErgoMax™** 24XSI

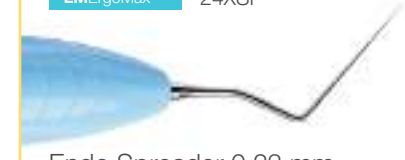

Endo Spreader 0.23 mm

- LMErgoSense®** 22ES
- LMErgoSense®** 22EST
- LMErgoMax™** 22XSI

ENDO PLUGGER

For the vertical plugging of root canal fillings.



Endo Plugger, 0.5–1.0 mm

- LMErgoSense®** 41-42ES
- LMErgoSense®** 41-42EST
- LMErgoMax™** 41-42XSI
- LMErgoNorm™** 41-42SI

Extract with confidence & luxate atraumatically

Familiar LM™ ergonomics and tactile sensitivity also in extracting instruments.

- Ergonomic design
 - comfortable to hold and easy to rotate
- Secure LM-ErgoTouch™ surface
 - light, non-slip grip
- Well-balanced
 - controlled grip
- Lightweight
 - less strain on the hands
- Clear color-coding
 - easy identification
- LM-DuraGradeMAX™ supersteel
 - strong and sharp blades
- Tight seam at blade juncture
 - hygienic

LM-LIFTOUT™



Luxating instruments for atraumatic extraction.

LM-LiftOut™ instruments allow tooth extraction to be done as atraumatically as possible. This is important in order to enable rapid healing and future implant placement. The thin tip of the instrument is introduced into the periodontal space and slowly advanced toward the apex of the root while twisting the instrument gently back and forth. LM-LiftOut™ is not suitable for use as an elevator.



STRAIGHT BLADE



2,5 mm

LM-LiftOut™ S2

812220

812220T



3 mm

LM-LiftOut™ S3

812230

812230T



4 mm

LM-LiftOut™ S4

812240

812240T



5 mm

LM-LiftOut™ S5

812250

812250T

CURVED BLADE



3 mm

LM-LiftOut™ C3

812430

812430T



4 mm

LM-LiftOut™ C4

812440

812440T



5 mm

LM-LiftOut™ C5

812450

812450T

LM-SLIMLIFT™

Luxating instruments designed for extremely narrow places.

LM-SlimLift™ instruments are especially designed to reach places you cannot go with traditional luxating instruments. It allows tooth extraction to be done as atraumatically as possible. This is important in order to enable rapid healing. The extremely fine tip of the instrument is introduced into the very narrow periodontal space. The instrument is slowly advanced towards the apex of the root while turning the instrument gently back and forth. The LM-SlimLift™ is not suitable for use as an elevator.



How to use

STRAIGHT BLADE**LM-SlimLift™ S3**

812630
 812630T

LM-SlimLift™ S4

812640
 812640T

LM-SlimLift™ S5

812650
 812650T

CURVED BLADE**LM-SlimLift™ C3**

812830
 812830T

LM-SlimLift™ C4

812840
 812840T

LM-SlimLift™ C5

812850
 812850T



The extremely slim design of blade allows for the most atraumatic extractions.

MAINTENANCE**Cassette**

810001

LM™ extraction instruments are supplied in a convenient cassette that protects both the instrument and the handler during the maintenance cycle. The cassette keeps the instrument from puncturing the sterilization pouch, so sterility is more likely to be assured.

- The cassette can be color-coded
- The cassette is compatible with LM-ServoMax™ tray

**Sharpening stone**

818002

Sharpen the extractors occasionally using the round sharpening stone as a simple and quick way to keep them in top condition and always ready to use.

- Suitable for LM-SlimLift™, LM-LiftOut™, LM-TwistOut™ and LM-RootOut™ instruments

LM-TWISTOUT™



Elevators for forceful extraction.

LM-TwistOut™ instruments are suitable for tooth extraction in situations where strong force or torque must be applied.



STRAIGHT
BLADE



3 mm

LM-TwistOut™ S3

814230
814230T



4 mm

LM-TwistOut™ S4

814240
814240T



5 mm

LM-TwistOut™ S5

814250
814250T

CURVED
BLADE



3 mm

LM-TwistOut™ C3

814430
814430T



4 mm

LM-TwistOut™ C4

814440
814440T



5 mm

LM-TwistOut™ C5

814450
814450T

ROOT TIP INSTRUMENTS



Root Tip Pick

- For the removal of root tip
- Suits well in all difficult to reach areas

LMErgoSense® 769-771ES
812210
LMErgoMax™ 769-771XSI
812210T

LM-RootOut™


- Designed for delicate root teasing procedures
- Not recommended for more substantial root tip removal
- Similar to Heidbrink root tip picks
- Ergonomically designed handle for a good grip



Mirror Handle, ortho

- Orthodontic mirror handle with Bausch-Verbiest design
- For measuring the sagittal overjet and vertical overbite
- Millimeter scale, total length 25 mm

LMErgoSense®	25-26ES
☞	25-26EST
LMErgoMax™	25-26XSI
LMErgoSense®	28-26ES
☞	28-26EST
LMErgoMax™	28-26XSI

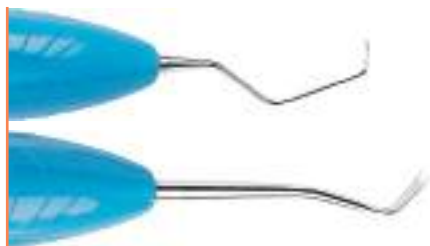

Ligature On-Off SL

- The explorer-like end (17CL) for both application and removal of ligatures. Ligatures can be lined up on the lower shank.
- The forked end (414) is particularly suitable for guiding a ligature ring over e.g. the fourth wing of a bracket.

LMErgoSense®	414-17CLES
☞	414-17CLEST
LMErgoMax™	414-17CLXSI



Plastic Ligature Applicator

- For placement of orthodontic rubber ligatures

LMErgoSense®	416-417ES
☞	416-417EST
LMErgoMax™	416-417XSI
LMErgoNorm™	416-417SI



Archwire Tucker

- For bending the end of an arch wire
- The flattened ball-ends have 1 mm holes, where the wire can easily be placed and turned.
- Long blade fits well even in tight places.

LMErgoSense®	412-413ES
☞	412-413EST
LMErgoMax™	412-413XSI



Ligature Tucker – Scaler U15

- For tucking of ligatures and arch wires (414)
- For removing of cement excesses following the removal of orthodontic appliance (150)


LMErgoSense®	414-150ES
☞	414-150EST
LMErgoMax™	414-150XSI



Ligature Director

- For tucking ligatures and arch wires and for placement of arch wire
- The forked end is particularly suitable for guiding a ligature ring over e.g. the fourth wing of a bracket
- The forked end tip is rotated 90°


LMErgoSense®	411-415ES
☞	411-415EST
LMErgoMax™	411-415XSI

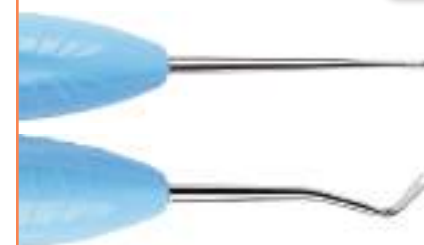

Ligature Tucker

- For tucking of ligatures and arch wires and for placement of arch wire
- The forked end is particularly suitable for guiding a ligature ring over e.g. the fourth wing of a bracket.

LMErgoSense®	414-415ES
☞	414-415EST
LMErgoMax™	414-415XSI
LMErgoNorm™	414-415SI

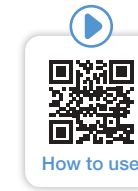



LM™ REPROCESSING RECOMMENDATIONS

To take best care of your premium quality products

DENTAL HAND INSTRUMENT MAINTENANCE WITH CARE

By following LM™ reprocessing recommendations, you can ensure your dental instruments to have both long life and high level performance. To take best care of your high-quality instruments pay attention to the recommendations and ensure that only validated procedures are used for cleaning and sterilization. For the hand instrument maintenance use LM-ServoMax™ tray system and LM-Servo™ cassettes in order to protect the instruments and avoid sharp injuries.

Always follow the regulations by local authorities.

LESS ADMINISTRATIVE CHORES AND MORE TIME FOR PATIENT CARE

LM Dental Tracking System™ automates traceability, logistics, and aids dental professionals in infection control documentation, streamlining processes and eliminating the necessity for labor-intensive manual tasks.

Implementation of LM Dental Tracking System™ enhances clinician and clinic manager knowledge that the essential reprocessing of materials and instruments has been done correctly, contributing to the safety of both clinicians and patients.



Scan the QR to read and download LM™ Reprocessing recommendations for hand instruments and accessories



LM Dental Tracking System™

FOR ENHANCE EFFICIENCY, SAFETY AND SAVINGS

LM™ INNOVATION

LM Dental Tracking System™ (LM DTS™) is developed to address everyday clinic challenges associated with dental material handling and traceability, including controlling clinic processes, instrument rotation, reuse, as well as inventory materials and consumables handling.

The unique system is created by combining scanning readers, proprietary server software and dental instruments and materials with built-in or external tags. Assets such as instruments, hand pieces, filling materials, implants, and more from various manufacturers can be wirelessly read in just a few seconds.

LM DTS™ automates traceability and dental asset handling processes, supports dental professionals in infection control documentation, and efficiently manages inventory materials and consumables handling. This streamlining of processes eliminates the need for labor-intensive manual tasks. The system allows you to turn vital information into action that makes a difference.

ENHANCED EFFICIENCY AND TIME SAVINGS

Automating instrument and material tracking, including documentation, significantly saves staff time previously spent on manual checks.

- Allocate additional time to patient care by minimizing administrative duties
- Assurance of instrument and material availability
- Paperless clinic

INFORMED DECISION-MAKING

Consistently improve practices by continually analyzing performance and monitoring real-time status.

- Ongoing performance improvement
- Instant dental inventory access
- Accessing cloud-based server from anywhere

ASSURED SAFETY CONFIDENCE

Enhances knowledge that the essential reprocessing of materials and instruments has been done correctly, contributing to the safety of both clinicians and patients.

- Automated checks and alerts
- Material traceability
- Maintenance and sterilization records
- Reprocessing procedure documentation



Scan the QR for more information at dentaltracking.com



LM-SERVO™ E INSTRUMENT CASSETTES



In LM-Servo™ E cassettes the instrument handles are clearly visible providing RFID tag readability and identification of individual instruments, though the tips are securely covered with safety shields preventing sharp injuries and protecting the blades.

The LM-Servo™ E cassette series is available with built-in RFID tagging. This option makes E cassettes excellent tool for improving traceability by enabling easy handling and tracking of instruments as sets. Below you find the product codes for RFID tagged cassettes marked with RFID symbol.

The cassettes can be coded with silicone code buttons (6000 – 6005) and attached to the LM-ServoMax™ with a special lock (6417).

LM-Servo™ 5E

- For five hand instruments
- Cassette is equipped with safety shields of different colors: grey, blue, green and yellow.
- Outer dimensions: 180 x 83 x 26 mm
- Inner dimensions: 174,8 x 78,6 x 21,4 mm




6551 grey
 6551T grey



6552 blue
 6552T blue



6553 green
 6553T green



6554 yellow
 6554T yellow

LM-Servo™ 8E

- For eight hand instruments
- Cassette is equipped with safety shields of different colors: grey, blue, green and yellow.
- Outer dimensions: 180 x 133 x 26 mm
- Inner dimensions: 174,8 x 128,9 x 21,4 mm




6581 grey
 6581T grey




6582 blue
 6582T blue



6583 green
 6583T green



6584 yellow
 6584T yellow

LM-SERVOMAX™ TRAY SYSTEM MD CE

LMServoMax™

Optimal work flow & rational maintenance.


LM-ServoMax™ tray system organises and rationalises the handling of dental instruments and accessories during the procedures and maintenance. It facilitates good infection control and saves time by reducing the handling of individual items. It also prolongs the life span of instruments and protects the personnel. Its versatile modules keep all instruments and accessories organised in a single tray. The system is adaptable to treatment specific and individual needs.



LM-ServoMax™ tray set
6410

- Includes: 1 set of universal trappers and 2 metal clips
- Dimensions: 288 x 185 x 29 mm

Rational maintenance with minimal handling.

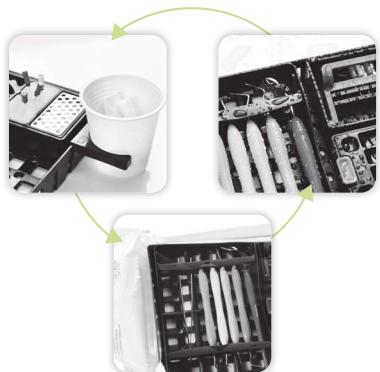

LM-ServoMax™ universal trapper,
2 pcs
6416



LM-ServoMax™ metal clips, 2 pcs
6415

- For locking LM-Servo™ 5 or LM-Servo™ 8 instrument cassette on the tray



LM-ServoMax™ standard case
6420



LM-ServoMax™ restoration / endo case
6430

LM-ServoMax™ foam, 20 pcs
6435

- For cleaning the endo files



LM-ServoMax™ ultrasonic tip case
6460

- For six ultrasonic tips



LM-ServoMax™ bur endofile case
6450

- For eight burs and fourteen files
- Includes adjustment and measuring scales for files



LM-ServoMax™ syringe holder
6495

- A stand attachable to the LM-ServoMax™ tray to place the anesthetic syringe safely in the needle guard






■ LM-ServoMax™ tray system modules have a universal design that complies the standard for dental tray dimensions. They are made of exceptionally durable special plastic which tolerate all common cleaning and sterilisation methods and agents (max 134 °C / 273 °F). Modules are sold empty.



LM-ServoMax™ waste cup holder
6490



LM-ServoMax™ instrument rack

-  6471 black
-  6472 orange
-  6473 blue
-  6474 red
-  6475 green

- Color-coded rails for attaching instruments to the LM-ServoMax™ tray
- Length adjustable according to the size of instruments
- Package 2 pcs / color

LM-SERVO™ CASSETTES



In LM-Servo™ cassettes hand instruments are locked securely but rotate freely for optimal cleaning and guaranteed infection control. The cassettes can be coded with silicone code buttons.

The cassettes can be attached to the LM-ServoMax™ tray with metal clips (6415).



Code Buttons

LM-Servo™ 5
6650

- For five hand instruments
- Dimensions: 180 x 86 x 28 mm



Code Buttons

LM-Servo™ 8
6680

- For eight hand instruments
- Dimensions: 180 x 136 x 28 mm

CODING PRODUCTS

Code buttons

-  6000 blue
-  6001 green
-  6002 orange
-  6003 red
-  6004 yellow
-  6005 brown

- For color-coding the LM-Servo™ and LM-Servo™ E cassettes
- Package 20 pcs / color



Code rings

-  5000 blue
-  5001 green
-  5002 orange
-  5003 red
-  5004 yellow
-  5005 black
-  5006 white
-  5007 brown
-  5008 lavender
-  5009 pink
-  5010 dark grey
-  5011 lime

- For color-coding hand instruments
- Package 50 pcs / color

Regular maintenance of dental hand instruments is the foundation of safe and high quality treatments.

In addition to cleanliness, instrument condition has a significant impact on the end result of a treatment as well as your and your patient's safety. When discussing hand instruments, a critical factor is the instrument blade. Only working with sharp instruments will ensure an optimal treatment result.

LM-RondoPlus™ is a simple and reliable sharpening machine, which is suitable for use with all hand instruments that can be resharpened. LM-RondoPlus™ sharpens all hand instruments easily, quickly and precisely.

Regular sharpening allows you to work with a sharp instrument, which reduces the amount of strain on your hand and enhances safety from one treatment to the next. The service life of a regularly sharpened instrument is also longer than one that is not maintained well.



Scan the QR for videos and more information.

LM™ INNOVATION

LMRondoPlus™ sharpening machine

360° rotating casing

Sharpening stone

LM-RondoPlus™ comes equipped with a ceramic 400 grit sharpening stone. A finer (600 grit) and a coarser (320 grit) stone are available as spare parts. The stone rotates in two directions. When the stone rotates away from the instrument, the sharpening effect is more delicate. For more sharpening power use a coarser stone that rotates toward the instrument. The direction of rotation is chosen from the rocker switch.

Instrument rest

The instrument rest has two positions. The lower position is for periodontal instruments and the upper position is for preparation instruments. The position is changed by lifting the rest slightly away from the sharpening stone and then up or down. The instrument rest must be as close as possible to the sharpening stone without actually touching it. Screw the control to adjust the distance between the rest and the stone.

Foot pedal connection

Use of the foot pedal is optional. When the foot pedal is connected, the sharpening stone will only rotate when you press down on the pedal. The rotational direction is still selected with the rocker switch.

Rocker switch
(direction control)

Testing sharpness

LM-RondoPlus™ is equipped with LM-Fingo™ sharpness test sticks, which can be mounted on either side of the machine or worn on the finger. Sharpness is easy to test by drawing the instrument lightly in its working position along LM-Fingo™. Sharp instrument will bite into the stick.

EASY SHARPENING

LM-RondoPlus™ is suitable for use with all hand instruments that can be resharpened. Thanks to its simple operation method, learning how to use the machine takes only a few moments.

The ceramic sharpening stone is mounted in a casing, which can rotate 360°. The instrument being sharpened is placed on the instrument rest next to the sharpening stone. The lower rest position is for periodontal instruments while the upper position is for preparation instruments. The hand holding the instrument remains still – the instrument is sharpened by rotating the casing.

With LM-RondoPlus™ the same basic technique is used, whether the instruments to be resharpened are periodontal or preparation instruments, or have convex or straight cutting edges.

LM-RondoPlus™ 8600

Package includes:

- Sharpening machine
- Foot pedal
- 400 grit sharpening stone
- 2 LM-Fingos™
- 5 emery tapes

Spare Parts

Grinding stone

- 8660 600 grit (white)
- 8640 400 grit, coarser (grey)
- 8632 320 grit, coarsest (black)

- 8630 Foot pedal
- 8550 Emery tapes, 5 pcs
- 8602 LM-Fingo™, 2 pcs
- 8606 LM-Fingo™, 6 pcs



On either side of the LM-RondoPlus™ there are LM-Fingo™ test sticks. Sharpness can be tested both with Fingo mounted on the LM-RondoPlus™ or your finger. The Fingo is made of an acrylic plastic and is sized to fit all users.

Easy to use

- Simple method - the blade is held in the correct sharpening position and the grinding stone rotates.
- Sharpening one instrument only takes a few seconds.

Versatile

- For both periodontal and preparation instruments.
- For both right or left handed.
- The direction of rotation can be selected.

Accurate and safe

- A dual-position instrument rest: The upper position is for preparation instruments and the lower for periodontal instruments.
- The original blade form is preserved.
- Low rotational speed: the blade will not heat up nor will its temper be affected.



LM-Activator™ - Early orthodontic treatment

LM-Activator™ is a prefabricated appliance for supporting healthy growth and development of the face and jaw. The ideal time for treatment with LM-Activator™ is in early mixed dentition when the first teeth are changing. By intervening early, the need for further treatment can be significantly reduced and even eliminated. Research indicates that malocclusions in deciduous dentition persist and tend to become more severe as the dentition develops. LM-Activator™ provides soft guidance to teeth and jaws when needed.

How does LM-Activator™ work?

- LM-Activator™ acts as a guide towards the correct position and place. Teeth and jaws are guided into correct positions and relation while teeth are changing in early mixed dentition.
- The forces of erupting teeth and the growth are utilized. The light and interceptive treatment supports the natural growth and development of teeth and jaws, and malocclusions are corrected.
- LM-Activator™ can be used at different ages. When the first milk tooth falls out, the opportunities for guiding both the teeth and jaws are wide.
- LM-Activator™ corrects and aligns sagittal and vertical relationships simultaneously - no separate phases are needed. The occlusion can be controlled in a 3-dimensional way.

Scan the QR and watch the video:
How does LM-Activator™ work



Gently guiding the growth of erupting teeth and growing jaws can be compared to trimming and supporting a young plant. Guiding to a healthy growth path early and along the growth is more gentle and the outcome will be more natural.



“The appliance has many advantages. Most importantly, early treatment with LM-Activator™ will stop the development of the malocclusions, thereby preventing the malocclusion to become more severe. For children, LM-Activator™ is easy to use. The treatment is easy and straightforward, and the treatment outcome more natural.”

Professor Emeritus Juha Varrelä, DDS, PhD, Specialist in Orthodontics, Institute of Dentistry, University of Turku, Finland

Benefits for A DENTIST

Corrects malocclusions effectively

Corrects typical malocclusions like excess Overjet, excess Overbite, Deep bite, Open bite, Class II and dental anterior Crossbite. Sagittal and vertical relationships can be corrected simultaneously.

✓ **Preventive.** Prevents the malocclusion from developing further.

✓ **Reduced need for treatment.** Can significantly reduce the need for further treatment.

✓ **Clinically and scientifically proven.** LM-Activator™ has been on the market since 2004 and the method is supported by evidence from cohort controlled and RCT research results.

Benefits for A YOUNG PATIENT

Comfortable to use

✓ **Made of medical-grade silicone.** Flexible by nature and comfortable for the patient.

✓ **Hygienic.** Daily rinse with water is enough. Can be immersed in boiling water.

✓ **Removable and used while sleeping.** Patient can eat any food normally, exercise and do any sport. No metal wire adjustments or emergency visits are needed due to broken metal parts.

✓ **Feels comfortable.** Anatomically designed and rounded slots.

✓ **Gentle on teeth.** Soft silicone is gentle on teeth and comfortable for the patient.

✓ **Simple and easy treatment.** No separate treatment phases but simultaneous alignment, leveling and anteroposterior correction.

Cost-effective treatment method

Reduced treatment complexity enabling increased savings and time saved for treating more patients.

✓ **Diverse.** Convenient assortment of models.

✓ **Cost-effective.** Appliance is affordable.

✓ **Prompt.** Treatment can be started immediately.

✓ **Efficient.** Light on chair-times as check-ups are short and not that frequent.

✓ **Collaborative.** The clinical work does not require full involvement by orthodontists.

Safe for young patients

✓ **No additives.** No coloring agents, phtalates, softeners, latex, bisphenol-A or fragrance agents.

✓ **Fully biocompatible.** The use is safe for the patient. The silicone material is tested in accordance with ISO 10993 requirement of 30 days permanent contact with mucous membranes.

✓ **No need for emergency visits** because of the soft silicone material.

✓ **Evidence based treatment method.** The treatment method is proven by research results, and the appliance has been used by a million of patients around the world.



INDICATIONS - when to use LM-Activator™?

LM-Activator™ has proven to be effective in treatment of different malocclusions at different dental stages. Evaluating the dentoalveolar, skeletal and functional characteristics is an essential part of selecting patients for LM-Activator™ treatment. A thorough case-by-case analysis of the patient's malocclusion and its severity as well as of the patients motivation is needed.

The ideal time for treatment with LM-Activator™ is in early mixed dentition when the first teeth are changing. Use of LM-Activator™ is recommended to align teeth, activate mandibular growth and expand the arch perimeter. The chart below provides a demonstrative guideline for case selection.

		Deciduous dentition	Early mixed dentition 5–8 years old	Late mixed dentition 9–13 years old	Early permanent dentition 13– years old
Excess overjet, excess overbite, deep bite	Severe		■	■	
	Moderate		■	■	■
	Mild		■	■	■
Anterior crowding (incisors and canines)	Severe		■	■	
	Moderate		■	■	■
	Mild		■	■	■
Dentoalveolar anterior crossbite of 1-2 teeth	Severe		■		
	Moderate		■		
	Mild		■		
Rotated anterior teeth	Severe		■		
	Moderate		■	■	
	Mild		■	■	■
Class II	Severe		■	■	
	Moderate		■	■	■
	Mild		■	■	■
Scissorsbite / Complete mandibular posterior lingual crossbite (i.e. mandibular posterior teeth completely on the lingual side of the maxillary posterior teeth)	Severe	■			
	Moderate	■			
	Mild	■			
Gummy Smile (Treatment start before the eruption of permanent maxillary central incisors)	Severe		■		
	Moderate		■		
	Mild		■		
Open bite	Severe		■	■	
	Moderate		■	■	■
	Mild		■	■	■

■ Use of LM-Activator™ recommended

■ Use of LM-Activator™ can be considered

Effect not proven for:

- Skeletal Class III
- Midline discrepancy > 3mm *)
- Very narrow upper arch *)
- Palatally impacted teeth *)
- Fully erupted anterior teeth that require torquing *)
- Teeth requiring intrusion

*) May be treated with combination treatment (LM-Activator™ together with other appliances such as quad helix)

Proven by research | References

[1] Keski-Nisula, K., Hernesniemi, R., Heiskanen, M., Keski-Nisula, L., & Varrelä, J. Orthodontic intervention in the early mixed dentition: A prospective, controlled study on the effects of the eruption guidance appliance. *AM J of Orthodontics and Dentofacial Orthopedics*, 2008;133(2), 254-260. <https://doi.org/10.1016/j.ajodo.2006.05.039>

[2] Keski-Nisula K., Keski-Nisula L., Varrelä J. Class II treatment in early mixed dentition with the eruption guidance appliance: effects and long-term stability. *European Journal of Orthodontics*, 2020;42(2):151-156. <https://doi.org/10.1093/ejo/cjz092>

[3] Keski-Nisula K., Keski-Nisula L., Salo H., Voipio K., Varrelä J. Dentofacial changes after Orthodontic Intervention with Eruption Guidance Appliance in the Early Mixed Dentition. *Angle Orthod* 2008;78(2):324-331. <https://doi.org/10.2319/012607-37.1>

[4] Myrland R., Dubland M., Keski-Nisula K., Kerosuo H. One year treatment effects of the eruption guidance appliance in 7- to 8-year -old children: a randomized clinical trial. *European Journal of Orthodontics*, 2015;37(2):128-134. <https://doi.org/10.1093/ejo/cju014>

[5] Myrland R., Keski-Nisula K., Kerosuo H. Stability of orthodontic treatment outcomes after 1-year treatment with the eruption guidance appliance in the early mixed dentition: A follow-up study. *Angle Orthod* 2019; 89(2):206-213. <https://doi.org/10.2319/041018-269.1>

[6] Kerosuo, H., Heikinheimo, K., Nyström, M., Väkiparta, M. Outcome and long-term stability of an early orthodontic treatment strategy in public health care. *European Journal of Orthodontics* 2013;35(2):183-189. <https://doi.org/10.1093/ejo/cjs087>

[7] Bishara S., B. Hoppens B., Jakobsen J., F. Kohout. Changes in the molar relationship between the deciduous and permanent dentitions: A longitudinal study. *AM J Orthodontics and Dentofacial Orthopedics* 1988;93(1):19-28. [https://doi.org/10.1016/0889-5406\(88\)90189-8](https://doi.org/10.1016/0889-5406(88)90189-8)

[8] Glasl B., Ludwig B., Schopf P. Prevalence and Development of KIG-relevant Symptoms in Primary School Students from Frankfurt am Main* Prävalenz und Entwicklung KIG-relevanter Befunde bei Grundschulern aus Frankfurt am Main*. *Journal of Orofacial Orthopedics / Fortschritte der Kieferorthopädie* 2006;67(6):414-423. <https://doi.org/10.1007/s00056-006-0615-8>

[9] Janson G. R., da Silva C. C., Bergersen E. O., Henriques J. F., Pinzan A. Eruption guidance appliance effects in the treatment of Class II, Division 1 malocclusions. *AM J of Orthodontics and Dentofacial Orthopedics* 2000;117(2), 119-129. [https://doi.org/10.1016/s0889-5406\(00\)70222-8](https://doi.org/10.1016/s0889-5406(00)70222-8)

The treating doctor is solely responsible and liable for diagnosis, treatment, and evaluating whether treatment with LM-Activator™, LM-Trainer™ or any other LM-Dental™'s product is appropriate for a specific patient, and for the outcome of any treatment with LM-Activator™, LM-Trainer™ or any other LM-Dental™'s product. Product used, and use of products shall be determined solely by doctor and evaluated for use with each individual patient where appropriate. Please refer to manufacturer's instructions for use.

TYPICAL MALOCCLUSIONS treated with LM-Activator™

Overjet ¹⁾



Before



After 1 year

Single tooth anterior crossbite ²⁾



Before



After 15 days

Deep bite with palatal impingement ¹⁾



Before



After 1 year

Crowding ³⁾



Before



After 1 year



“Our results showed distinct improvements in overjet, overbite, sagittal molar relationship, and crowding in the treated subjects.”

Myrlund R., et al. One year treatment effects of the eruption guidance appliance in 7- to 8-year-old children: a randomized clinical trial. *European Journal of Orthodontics*, 2015;37(2):128-134. <https://doi.org/10.1093/ejo/cju014>



“An advantage of the eruption guidance appliance is that it not only guides the eruption of the teeth but also simultaneously acts on the transversal, sagittal, and vertical relationships of both dental arches.”

Keski-Nisula, K., et al. J. Orthodontic intervention in the early mixed dentition: A prospective, controlled study on the effects of the eruption guidance appliance. *AM J of Orthodontics and Dentofacial Orthopedics*, 2008;133(2), 254-260. <https://doi.org/10.1016/j.ajodo.2006.05.039>



“As a result of the EGA treatment, the sagittal relationship was corrected from Class II to Class I in 86% of the cases during the active treatment, and it showed further improvement post-treatment. At the age of 16.7 years, 98% of the treatment children, who all had a Class II relationship at the onset of the early mixed dentition, had a Class I relationship.”

Keski-Nisula K., et al. Class II treatment in early mixed dentition with the eruption guidance appliance: effects and long-term stability. *European Journal of Orthodontics*, 2020;42(2):151-156. <https://doi.org/10.1093/ejo/cjz092>

Clinical pictures courtesy of ¹⁾ and ²⁾ Dr. Rita Myrlund and Dr. Mari Dubland, University of Tromsø, Norway. ³⁾ Dr. Gioacchino Pellegrino, Management of eruption problems with preformed silicone orthodontic devices. Poster presented at 21. Jahrestagung der Deutsche Gesellschaft für Kinderzahnheilkunde 25.09 - 27.09 2014; Freiburg, Germany

LM-Activator™ MODEL ASSORTMENT features

With LM-Activator™ alignment, leveling and anteroposterior correction can be done at the same time without separate phases. A large assortment of models makes it possible to select the appropriate appliance for the patient. Treatment complexity for both the clinician and the patient can be reduced.

The two basic LM-Activator™ model features are: **LOW** and **HIGH**. These are available with shorter or longer molar sections: **SHORT** and **LONG**. The appliances are combinations of these features. LM-Activator™ 2 models available also for two arch widths: **NARROW** or **WIDE**. **REINFORCED** model with hard insert in the incisal area is available. We recommend it for patients over 8 years of age only. The choice of model should always be based on an examination and diagnosis done by or under supervision of a specialist in orthodontic treatment.



LOW ●●

Model with basic appliance thickness is applicable to many different cases.

HIGH ●●

Model with thicker premolar and molar area. Specifically designed for treating skeletal and dentoalveolar open bite cases.



SHORT ●●

Model with a shorter molar section for patients whose second molars have not yet erupted.

LONG ●●

Model with a longer molar section for patients whose second molars have erupted.



NARROW


WIDE

There are two different arch widths of all basic models. The most commonly used arch width is **NARROW** with standard or regular arch. The **WIDE** model has a wider and shorter arch, and is especially designed for Asian markets.



REINFORCED

LM-Activator™ 2 models are available with incisal reinforcement. This model has a rigid reinforcement of harder material in the incisal area and is designed especially for skeletal deep bite cases. The model is also suitable for other patients who benefit from a firm surface and improved durability. Reinforced model is recommended for patients over 8 years old.

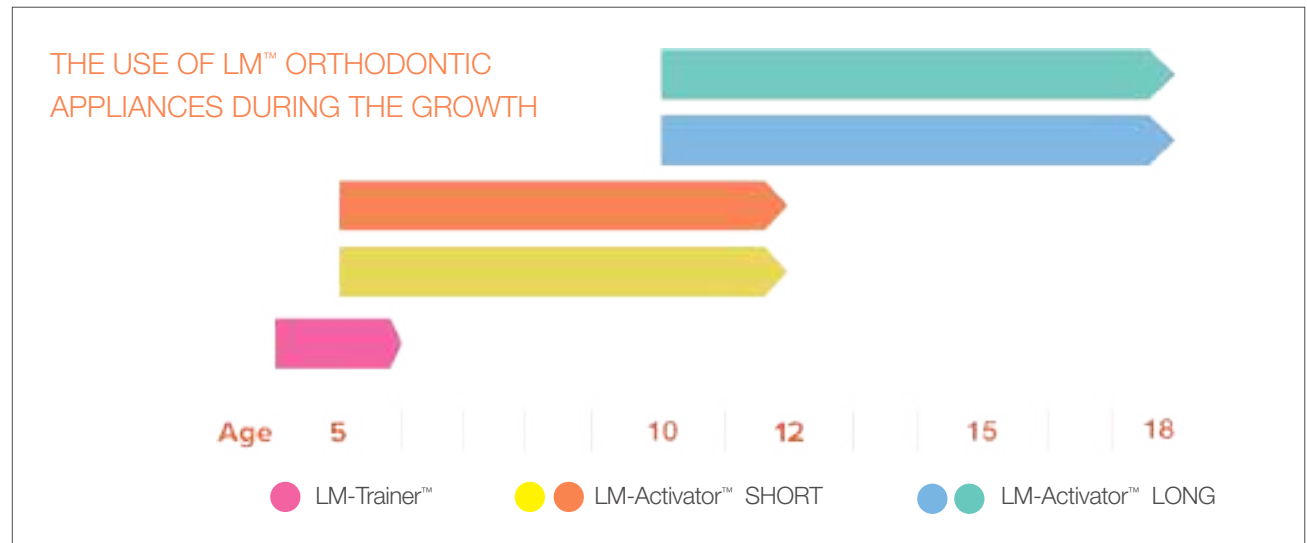

Selecting THE SIZE with LM-OrthoSizer™

LM-OrthoSizer™ is a tool for aiding the selection of the appropriate LM-Activator™ size. The size is indicative and LM-Activator™ should always be fitted to ensure the correct size. There are two models: one for LM-Activator™ (9400) and another for LM-Activator™ 2 (9402). Both are compatible with LM™ mirror handles (25, 28).



LM-Activator™ 2 OrthoSizer™ (9402) is used with LM-Activator™ 2 models. Separate measures for mandible and maxilla improve the accuracy.


LM-OrthoSizer™ (9400) is used with LM-Activator™ models.

LM-Activator™ 2 OrthoSizer™ measures the distance across the incisors from the distal surface of the left lateral incisor to the right lateral incisor.

Place the elevated marker between the left maxillary lateral incisor and canine. Ensure that the text *Maxilla* is at the upper jaw and the text *Mandibula* at the lower jaw.

Read the scale between the right maxillary lateral incisor and canine (i.e. at the mesial surface of the canine). Sizes 10, 20, 30... are indicated with larger markers with a dot whereas sizes 15, 25, 35... are indicated with smaller markers. (In the right side photo the size indication is 55.) If you wish to get an indication of the size separately for the mandible, use the LM-Activator™ 2 OrthoSizer™ in the same way in the mandible.

In case of crowding, consider a larger appliance and in case of diastemas, consider a smaller appliance. Fit LM-Activator™ on the patient and double check that the teeth sit properly in the slots both in the maxilla and in the mandible.

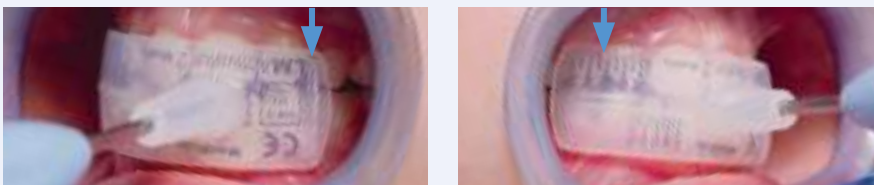

Fitting

Fitting the appliance to the patient is the most important step of size selection.



1. No crowding and no risk of crowding

Choose the size that matches the dentition. The canine rests at the bottom of the slot and no distalizing force is placed on the canine. If there is no crowding or risk of crowding, the size in the photo is correct. If there is crowding or crowding is expected, consider a larger size to enable expansion of the arch perimeter. In case of diastema consider choosing a smaller size.



2. Crowding or risk of crowding

This size is larger than in Figure 1.

The canine is guided by LM-Activator™ towards the bottom of its' slot and applies a force that pushes the canine distally. This size is appropriate if crowding is present or expected and arch perimeter expansion is needed.



3. The size is too large

The appliance is too large. The canine is guided against a ridge between two slots and the appliance does not guide teeth properly. Choose a smaller size.



At follow-up check-ups with the patient, ensure that the size is still correct.

Re-evaluate the size of the appliance when the permanent maxillary and mandibular lateral incisors begin to erupt. It may be necessary to switch to a larger appliance.



Scan the QR and watch the video:
Selecting the size and fit of LM-Activator™

General steps on HOW TO CHOOSE the LM-Activator™ MODEL

The selection chart below provides some general basic guidelines on choosing the right LM-Activator™ appliance for the patient. It is solely the responsibility of the treating doctor to evaluate and determine which model is best suitable for the individual patient.



Scan the QR and watch video: Guidelines on how to choose the LM-Activator™ model

STEP 1

Malocclusion type?

- Increased overjet
- Increased overbite/deep bite
- Anterior crowding
- Dentoalveolar anterior crossbite
- Rotated anterior teeth
- Class II
- Scissorsbite
- Complete mandibular posterior lingual crossbite
- Gummy smile



Anterior open bite
(Use High model of LM-Activator™)

Growth pattern / Facial type?



- Low angle Short face Hypodivergent Brachyfacial
- Neutral angle Neutral growth pattern Mesofacial
- High angle Long face Hyperdivergent Dolichofacial

STEP 2

Have the permanent second molars erupted?



STEP 3

Dental arch width?

- Narrow / Regular → NARROW
- Wide → WIDE

STEP 4

Is reinforcement needed in incisal area?

- Yes → REINFORCED

“We can control the occlusion in a three dimensional way with LM-Activator™. We are able to control overbite, overjet, the lateral contacts of the lower and upper arch and to modify the sagittal growth.”



Dr. Gioacchino Pellegrino, DDS, Specialist in Orthodontics, Italy

The treating doctor is solely responsible and liable for diagnosis, treatment, and evaluating whether treatment with LM-Activator™, LM-Trainer™ or any other LM-Dental™'s product is appropriate for a specific patient, and for the outcome of any treatment with LM-Activator™, LM-Trainer™ or any other LM-Dental™ product. Product used, and use of products shall be determined solely by doctor and evaluated for use with each individual patient where appropriate. Please refer to manufacturer's instructions for use.

The LM-Activator™ model features are: LOW and HIGH. These are available with shorter or longer molar sections: SHORT and LONG. The appliances are combinations of these features. The LM-Activator™ 2 models are available with incisal reinforcement for increased durability, REINFORCED. We recommend it for patients over 8 years of age only.

LM-Activator™ 2 models have two alternative arch widths. The most commonly used arch width is NARROW with standard or regular arch. The WIDE model has a wider and shorter arch, and is especially designed for Asian markets. It is suitable also for some Latin American markets. The choice of model should always be based on an examination and diagnosis done by or under supervision of a specialist in orthodontic treatment.



LM-Activator™ 2 product range MD CE

LM-Activator™ 2 Low Short ●

Model with a shorter molar section when the second molars have not yet erupted. Delivered in a yellow box.

SIZE	PRODUCT CODES			
	NARROW	REINFORCED	WIDE	REINFORCED
35	94235LSN	94235LSNR	94235LSW	94235LSWR
40	94240LSN	94240LSNR	94240LSW	94240LSWR
45	94245LSN	94245LSNR	94245LSW	94245LSWR
50	94250LSN	94250LSNR	94250LSW	94250LSWR
55	94255LSN	94255LSNR	94255LSW	94255LSWR
60	94260LSN	94260LSNR	94260LSW	94260LSWR
65	94265LSN	94265LSNR	94265LSW	94265LSWR
70	94270LSN	94270LSNR	94270LSW	94270LSWR

LM-Activator™ 2 Low Long ●

Model with a longer molar section when the second molars have erupted. Delivered in a blue box.

SIZE	PRODUCT CODES			
	NARROW	REINFORCED	WIDE	REINFORCED
35	94235LLN	94235LLNR	94235LLW	94235LLWR
40	94240LLN	94240LLNR	94240LLW	94240LLWR
45	94245LLN	94245LLNR	94245LLW	94245LLWR
50	94250LLN	94250LLNR	94250LLW	94250LLWR
55	94255LLN	94255LLNR	94255LLW	94255LLWR
60	94260LLN	94260LLNR	94260LLW	94260LLWR
65	94265LLN	94265LLNR	94265LLW	94265LLWR
70	94270LLN	94270LLNR	94270LLW	94270LLWR

LM-Activator™ 2 High Short ●

Model with a shorter molar section. The High model is thicker in the region of second premolars and molars. It is designed for treating skeletal and dentoalveolar Open bite cases. Delivered in an orange box.

SIZE	PRODUCT CODES	
	NARROW	WIDE
35	94235HSN	94235HSW
40	94240HSN	94240HSW
45	94245HSN	94245HSW
50	94250HSN	94250HSW
55	94255HSN	94255HSW
60	94260HSN	94260HSW
65	94265HSN	94265HSW
70	94270HSN	94270HSW

LM-Activator™ 2 High Long ●

Model with a longer molar section. The High model is thicker in the region of second premolars and molars. It is designed for treating skeletal and dentoalveolar Open bite cases. Delivered in a green box.

SIZE	PRODUCT CODES	
	NARROW	WIDE
35	94235HLN	94235HLW
40	94240HLN	94240HLW
45	94245HLN	94245HLW
50	94250HLN	94250HLW
55	94255HLN	94255HLW
60	94260HLN	94260HLW
65	94265HLN	94265HLW
70	94270HLN	94270HLW

LM-Activator™ product range MD CE

LM-Activator™ Low Short ●

Model with a shorter molar section when the molars have not yet erupted. Delivered in a yellow box.

SIZE	NARROW
10	94010LS
15	94015LS
20	94020LS
25	94025LS
30	94030LS
35	94035LS
40	94040LS
45	94045LS
50	94050LS
55	94055LS
60	94060LS
65	94065LS
70	94070LS



LM-Activator™ Low Long ●

Model with a longer molar section when the second molars have erupted. Delivered in a blue box.

SIZE	NARROW
40	94040L
45	94045L
50	94050L
55	94055L
60	94060L
65	94065L
70	94070L



LM-Activator™ High Short ●

Model with a shorter molar section. The High model is thicker in the region of second premolars and molars. It is designed for treating skeletal and dentoalveolar Open bite cases. Delivered in an orange box.

SIZE	NARROW
20	94020HS
25	94025HS
30	94030HS
35	94035HS
40	94040HS
45	94045HS
50	94050HS
55	94055HS
60	94060HS
65	94065HS
70	94070HS



LM-Activator™ High Long ●

Model with a longer molar section. The High model is thicker in the region of second premolars and molars. It is designed for treating skeletal and dentoalveolar Open bite cases. Delivered in a green box.

SIZE	NARROW
45	94045H
50	94050H
55	94055H
60	94060H
65	94065H
70	94070H



LM-Trainer™ products MD CE

LM-Trainer™ can be used in deciduous dentition, e.g. before LM-Activator™ treatment. It can also be used in functional education and for habit correction, e.g. reverse swallowing and mouth breathing, when harmful oral habits risk causing malocclusions.



LM-Trainer™ Small

94100S

- Slots for incisors, lip bumper
- Smaller and softer than LM-Trainer™ Medium



LM-Trainer™ Medium

94100T

- Slots for incisors and cuspids
- Lip bumper



LM-Trainer™ 2 Medium

94100T2

- Bigger than LM-Trainer™ Medium
- High gloss surface
- Slots for incisors and cuspids
- High labial shields and shaped lingual flanges
- Enhanced lip bumper



LM-Trainer™ Braces

94100TB

- Used in conjunction with brackets and wires
- Brings the mandible forward for Class II correction
- Prevents irritation and wounding of soft tissue
- Lip bumper

LM-Activator™ motivational materials

Patient compliance is essential for successful treatment. The first month is crucial in adopting the habit to use the appliance, because creating a habit in general requires 3-4 weeks of regularity. To motivate the patient and his/her family, age appropriate material is available for both preadolescents and adolescents. It is important to make sure the patient and the family understands that irregular use may be equivalent with no treatment at all.



LM-Activator™ patient motivation path with stickers M1063A

- The language-independent motivation path is used to motivate the daily use of the LM-Activator™ appliance. A young patient is supposed to place a sticker on the calendar path every time she/he wears the appliance (day or night).



LM-Activator™ Clinic brochure "Growing happy smiles" 1061EN

- Take-away brochure for clinics to introduce and promote LM-Activator™ treatment for parents. English.



How to use LM-Activator™ Patient Guidance brochure M1062EN

- Cartoon styled guidance for parents and patients on how to use LM-Activator™.
- Includes also a calendar for tracking daytime and night-time use. English.



LM-Activator™ Clinic poster M1059EN

- For promotion Early orthodontic treatment with LM-Activator™ in a waiting room at a dental clinic. English.



For more information on early orthodontics and LM-Activator™ visit www.lm-activator.com





Making dental care safe, efficient and pleasant.



LM Dental™

 LM-Instruments Oy

Norrbyn rantatie 8 | FI-21600 Parainen | Finland

Tel. +358 2 4546 400

info@lm-dental.com | www.lm-dental.com

Making dental care safe, efficient and pleasant.

